


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

1-(1-Naphthyl)-ethanol

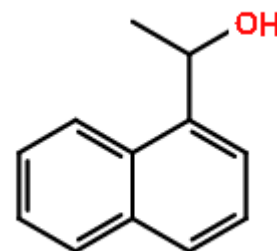

CHROMATOGRAM

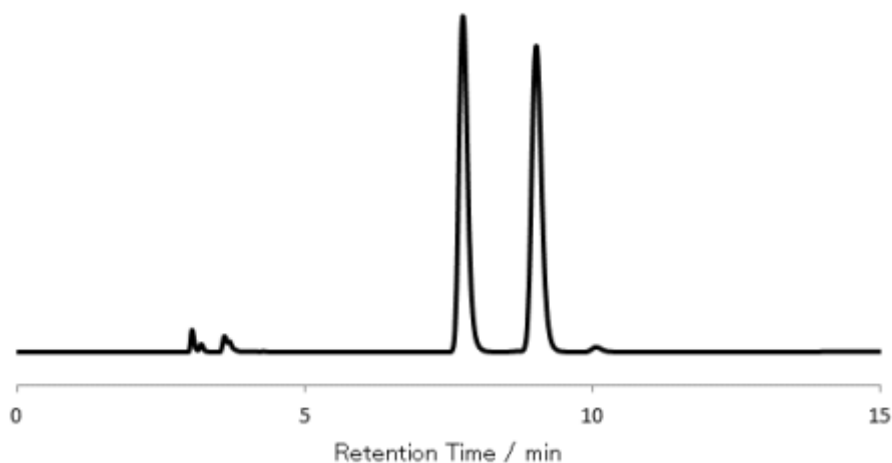

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IK		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol = 90 / 10		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	5	Sample Concentration	0.5 mg/mL in Mobile phase
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.7	1.55	1 - 2	4.11
2	9.0	1.97		