


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Omeprazole

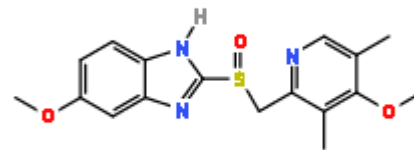

CHROMATOGRAM

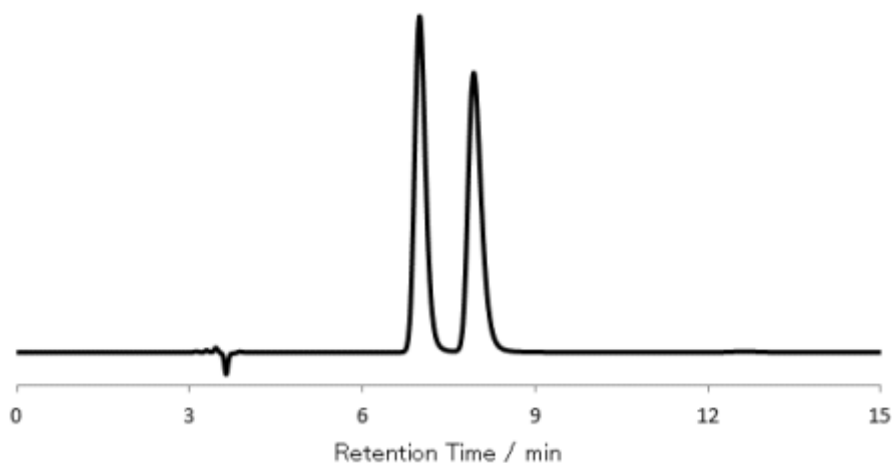

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IK		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 60 / 40 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	5	Sample Concentration	1 mg/mL in Mobile phase
Detection	UV-VIS 302 nm nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.4	1.46	1 - 2	6.93
2	9.7	2.22		
3				