


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Cetirizine

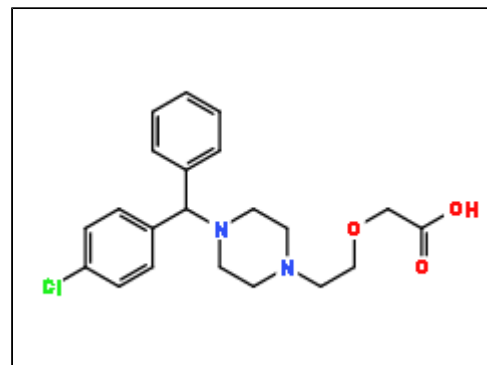

## CHROMATOGRAM

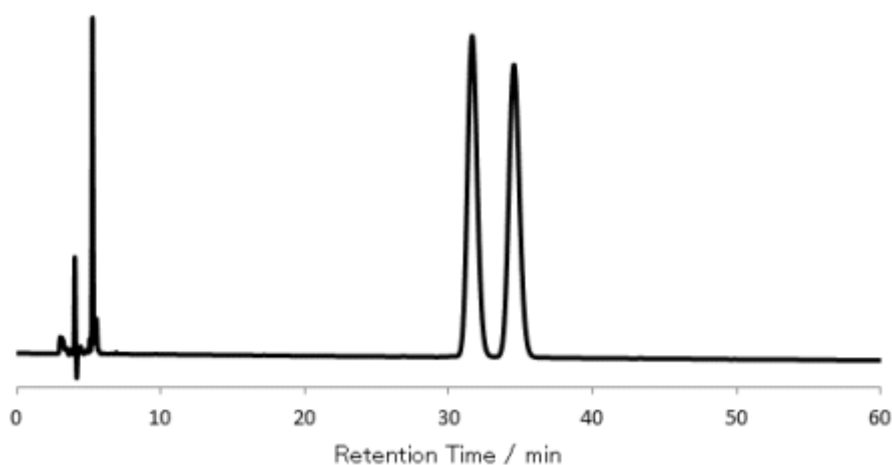

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IK-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / acetic acid / diethylamine = 70 / 30 / 0.1 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	5	Sample Concentration	1 mg/mL in Mobile phase
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	31.6	9.54	1 - 2	3.19
2	34.5	10.51		
3				