


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Bifonazole

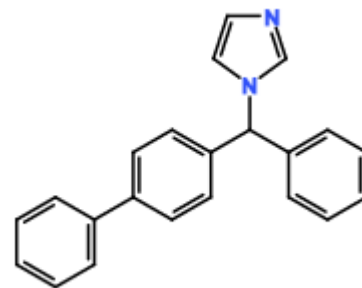

## CHROMATOGRAM

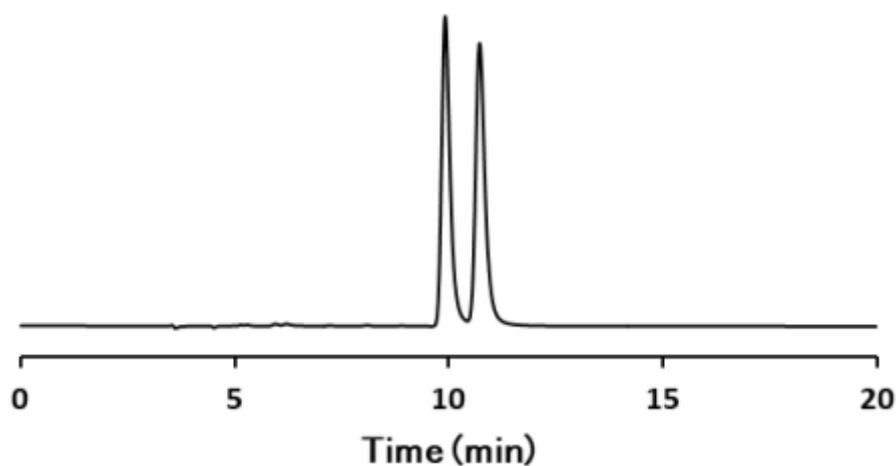

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IK		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	acetonitrile / diethylamine = 100 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	5	Sample Concentration	1.0 mg/mL
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.9	2.30	1 - 2	2.10
2	10.7	2.57		