


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

4-Benzoyloxy-2-azetidinone

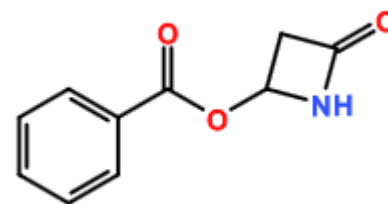

## CHROMATOGRAM

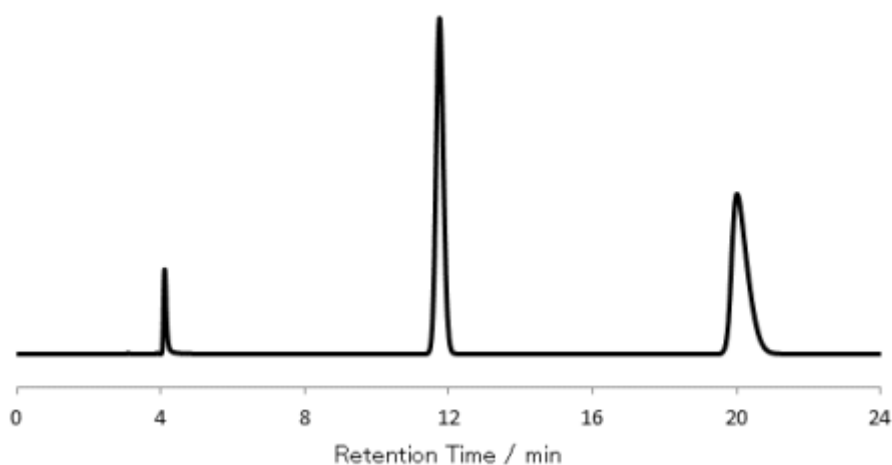

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IN		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol = 80 / 20		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	5	Sample Concentration	1.0 mg/mL in Mobile phase
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	11.7	2.91	1 - 2	13.48
2	20.0	5.66		