


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Clorprenaline

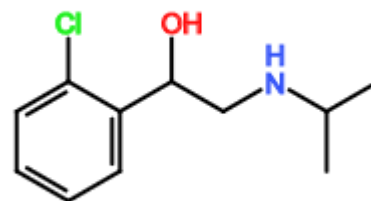

CHROMATOGRAM

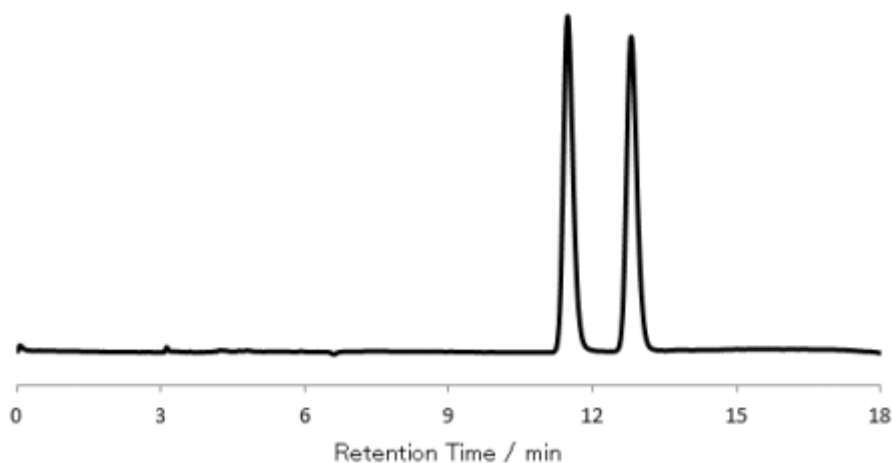

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IN		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 100 / 1 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	5	Sample Concentration	1.0 mg/mL in Mobile phase
Detection	UV-VIS 262 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	11.5	2.82	1 - 2 1.16	3.32
2	12.8	3.26		