


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Propafenone

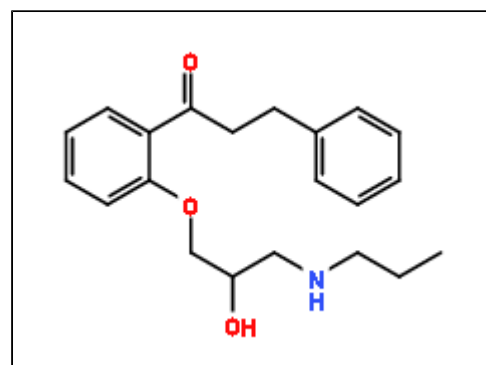

CHROMATOGRAM

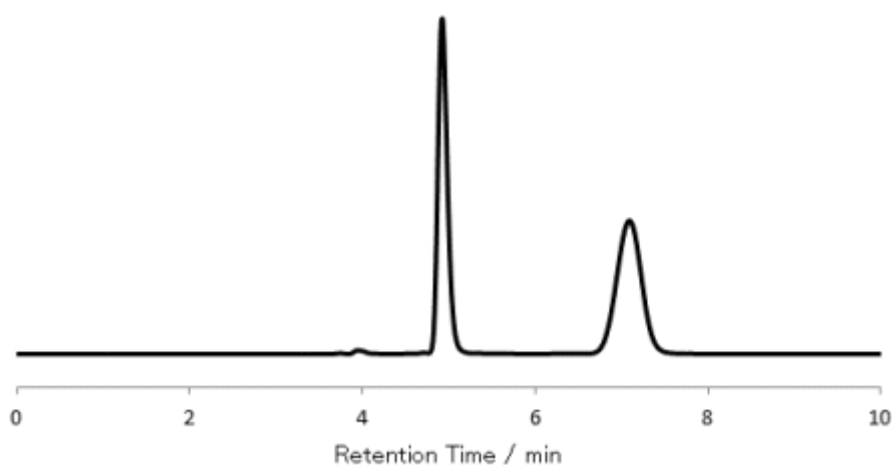

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IN		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	acetonitrile / ethanolamine = 100 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	10	Sample Concentration	1.0 mg/mL in Mobile phase
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.9	0.64	1 - 2	2.13
2	7.1	1.36		5.71