


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

4-Benzoyloxy-2-azetidinone

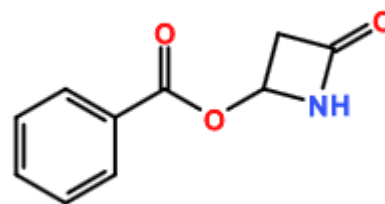

CHROMATOGRAM

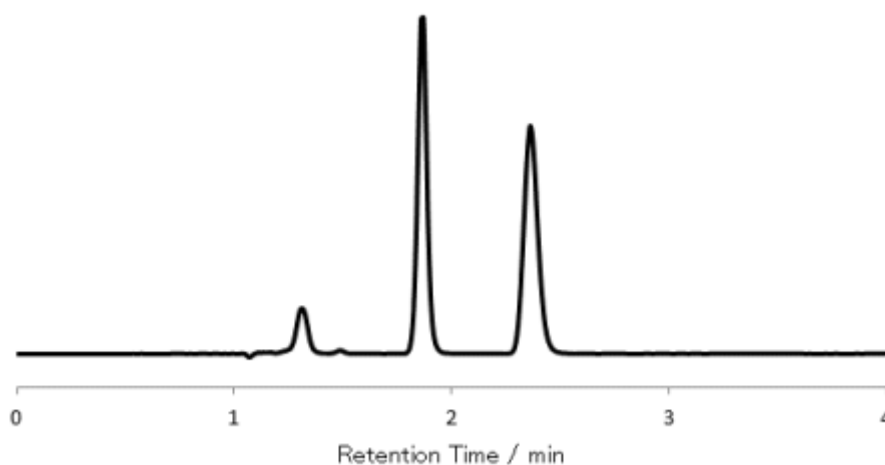

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IN		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	CO ₂ / methanol / diethylamine = 70 / 30 / 0.03		
Chromatographic Mode	SFC		
Flow rate (mL/min)	3.0	Temperature (°C)	40
Injection Volume (µL)	5	Sample Concentration	1.0 mg/mL in Mobile phase
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	1.9	0.68	1 - 2	5.10
2	2.4	1.13		