


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

(±)-N,N'-Bis(alpha-methylbenzyl) sulfamide

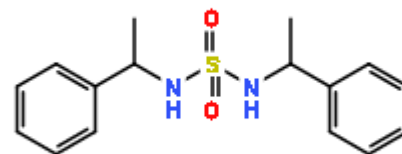

CHROMATOGRAM

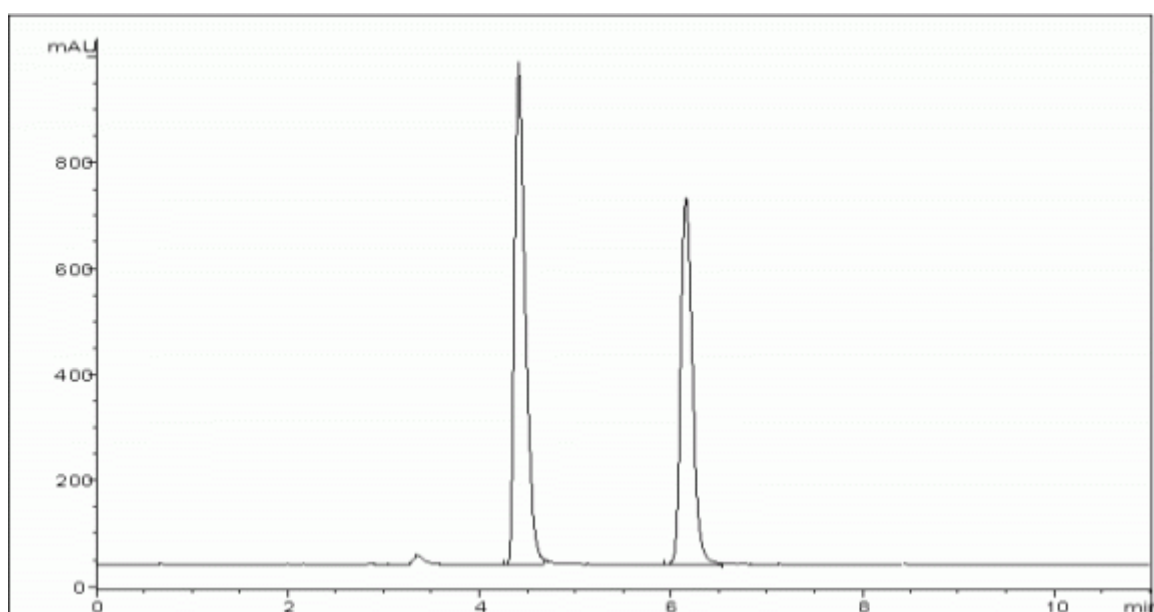

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	dichloromethane = 100		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.4	0.44	1 - 2	8.27
2	6.2	1.01		