


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

(4R,5S/4S,5R)-1,5-Dimethyl-4-phenyl-2-imidazolidinone

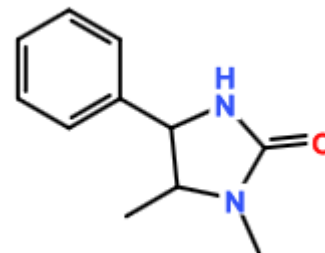

CHROMATOGRAM

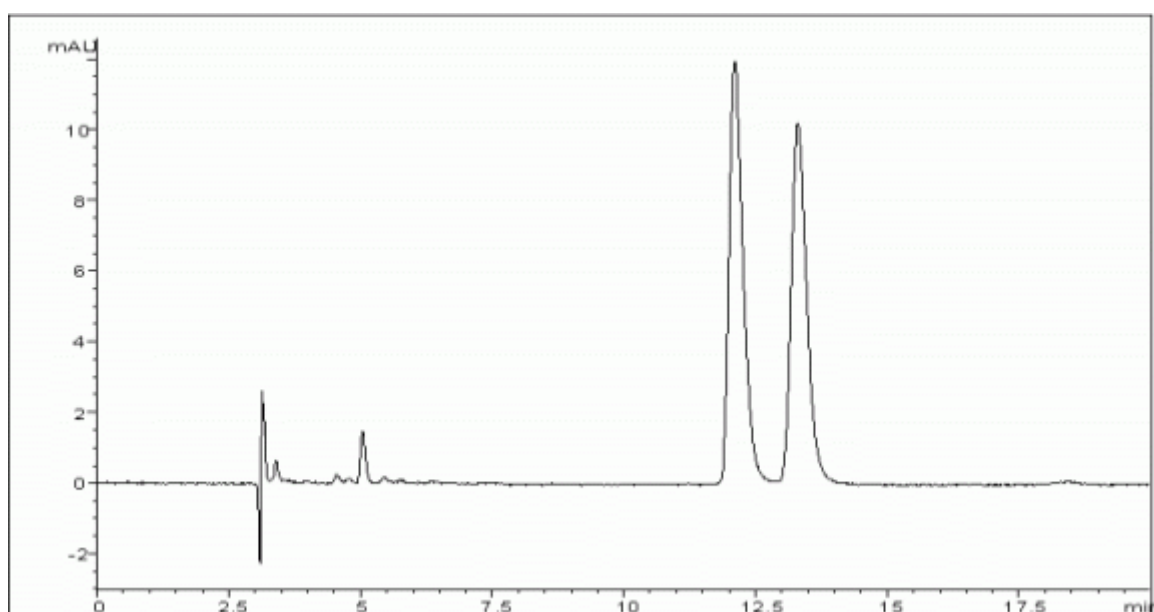

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethyl acetate = 20 / 80		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	12.1	2.95	1 - 2	2.40
2	13.3	3.33		