


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

2-Bromo-1-Indanone

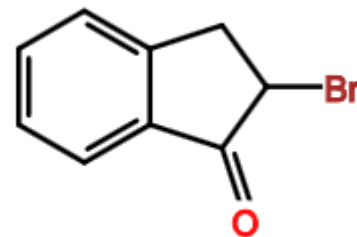

## CHROMATOGRAM

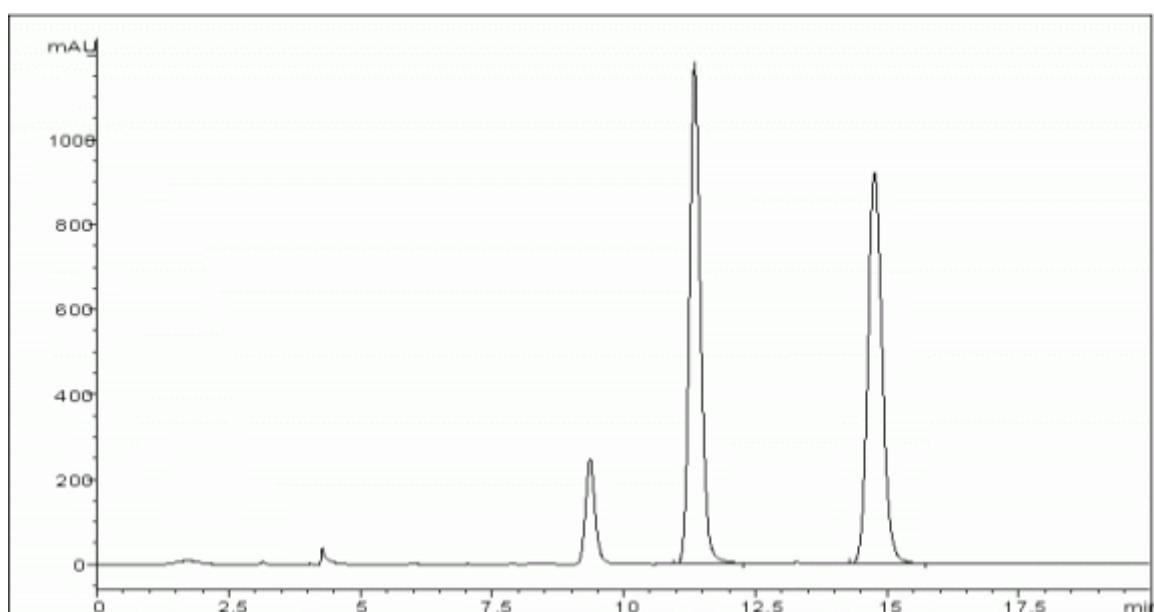

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol = 85 / 15		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	11.3	2.74	1 - 2	8.04
2	14.8	3.86		