


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

9,10-Dihydro-2-methyl-4H-benzo-[5,6] cyclohept[1,2-d]oxazol-4-ol

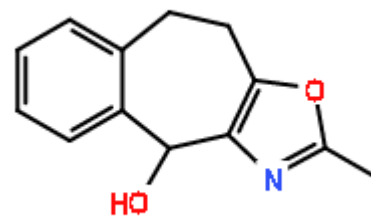

## CHROMATOGRAM

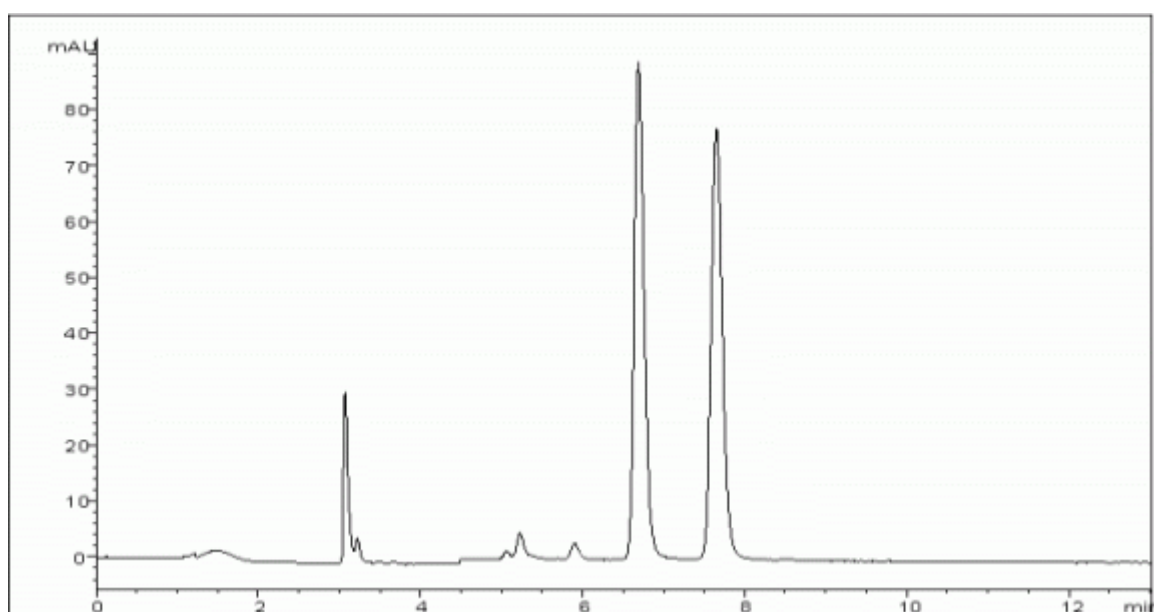

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethyl acetate / ethanolamine = 60 / 40 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.7	1.18	1 - 2	4.13
2	7.7	1.49		