


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

beta-Alanyl-DL-leucine

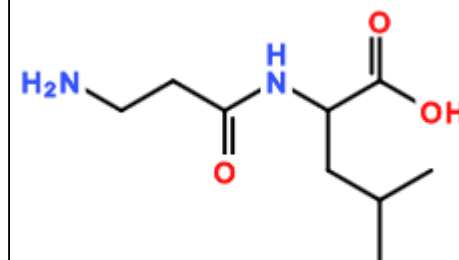

## CHROMATOGRAM

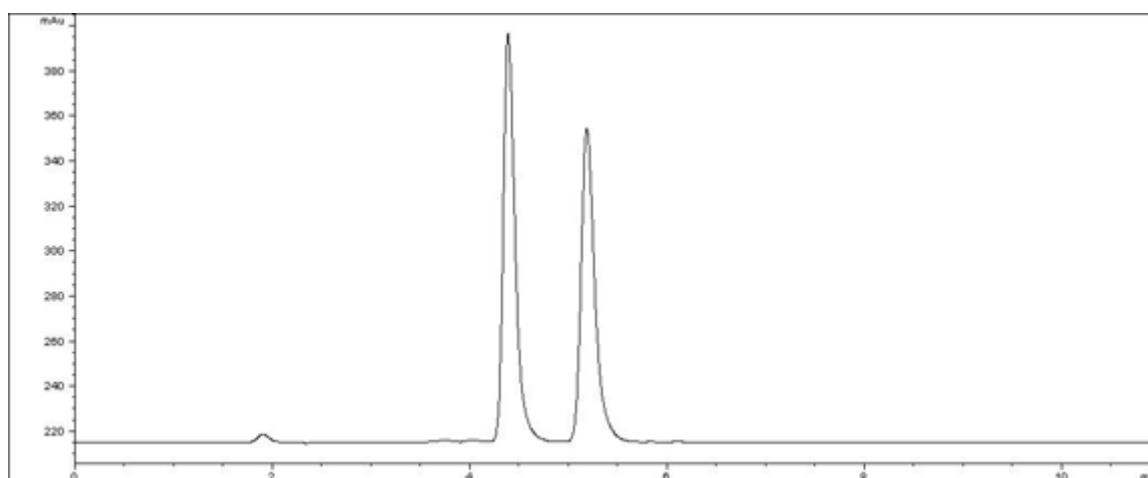

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	3 x 150 mm (3 µm)		
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / water = 98 / 2		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.0015 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.4	1.43	1 - 2	3.59
2	5.2	1.87		