

## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Cisatracurium (racemic)



## CHROMATOGRAM

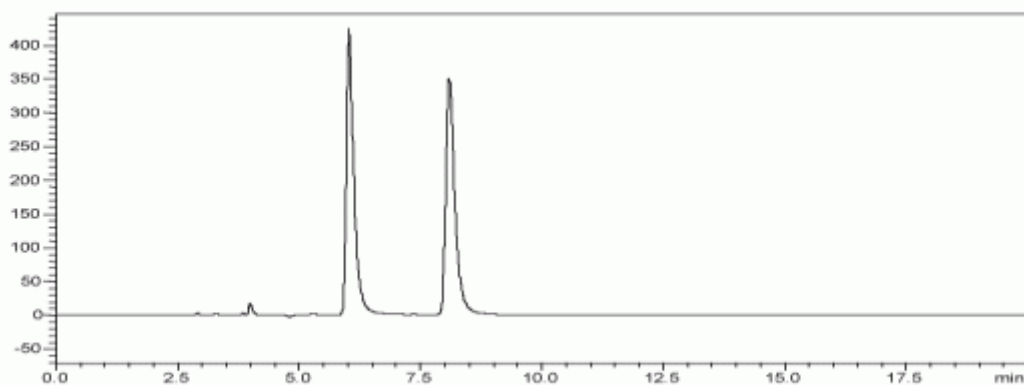

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / methanol / diethylamine = 80 / 10 / 10 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 240 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.0	1.00	1 - 2	7.04
2	8.1	1.70		