


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Methaqualone

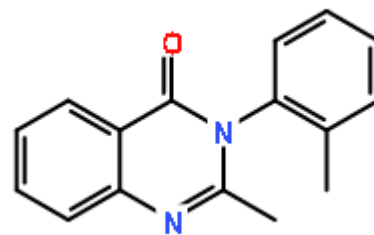

## CHROMATOGRAM

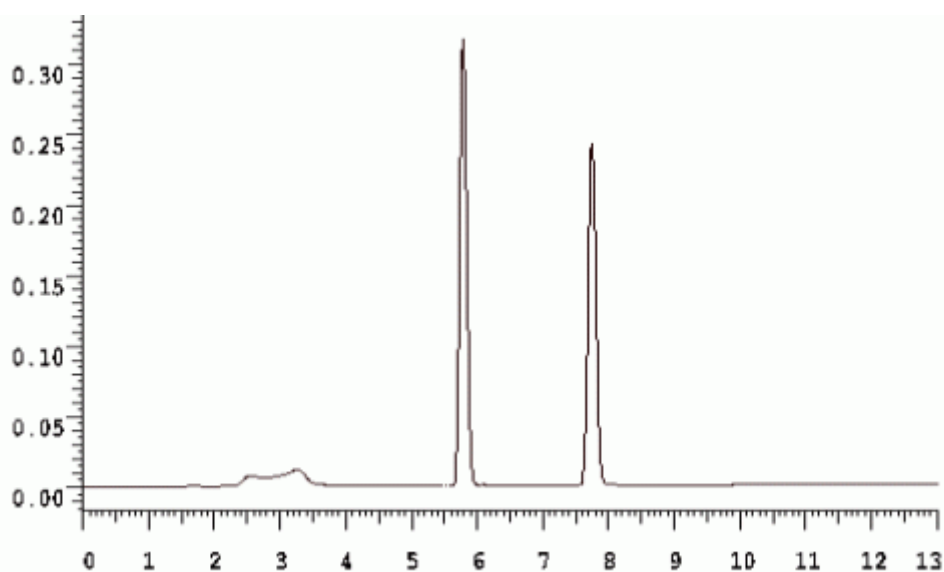

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	toluene / n-hexane / ethanol = 70 / 25 / 5		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	0.8	Temperature (°C)	25
Injection Volume (µL)	20 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.8	0.54	1 - 2	9.28
2	7.7	1.06		