


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Methyl-(phenylsulfinyl)-acetate

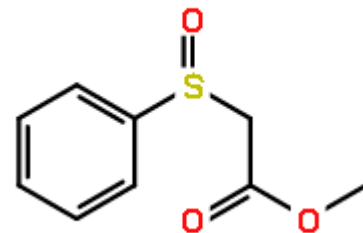

CHROMATOGRAM

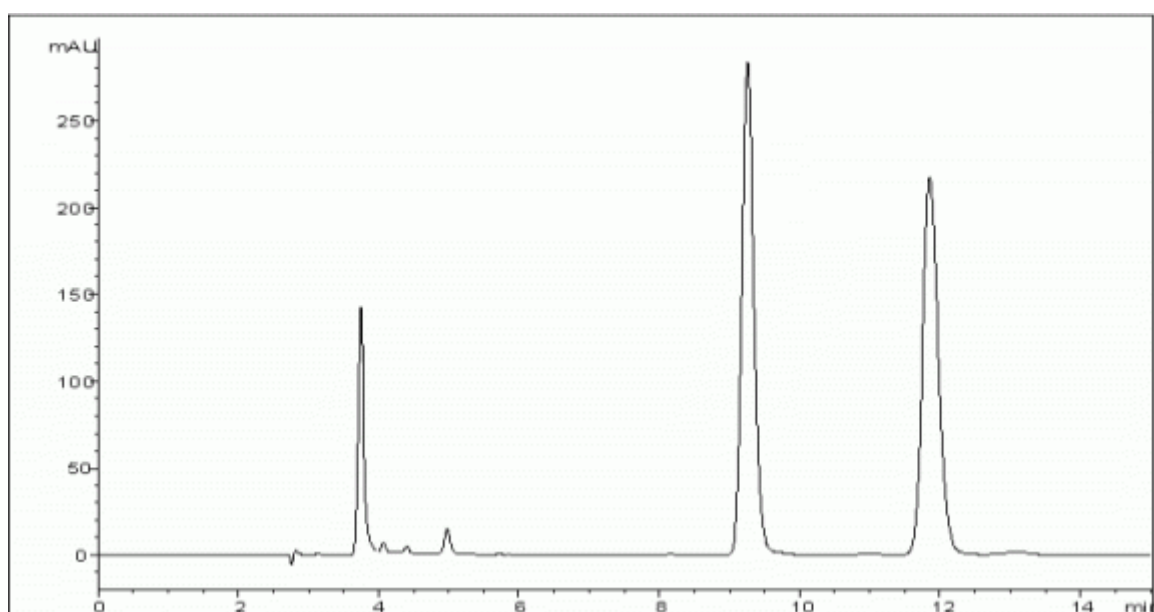

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-heptane / ethanol = 70 / 30		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 250 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.2	2.08	1 - 2	7.26
2	11.9	2.95		