


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Oxprenolol

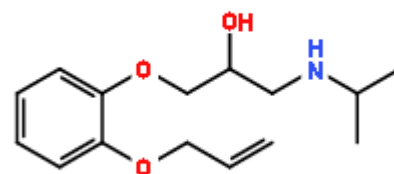

CHROMATOGRAM

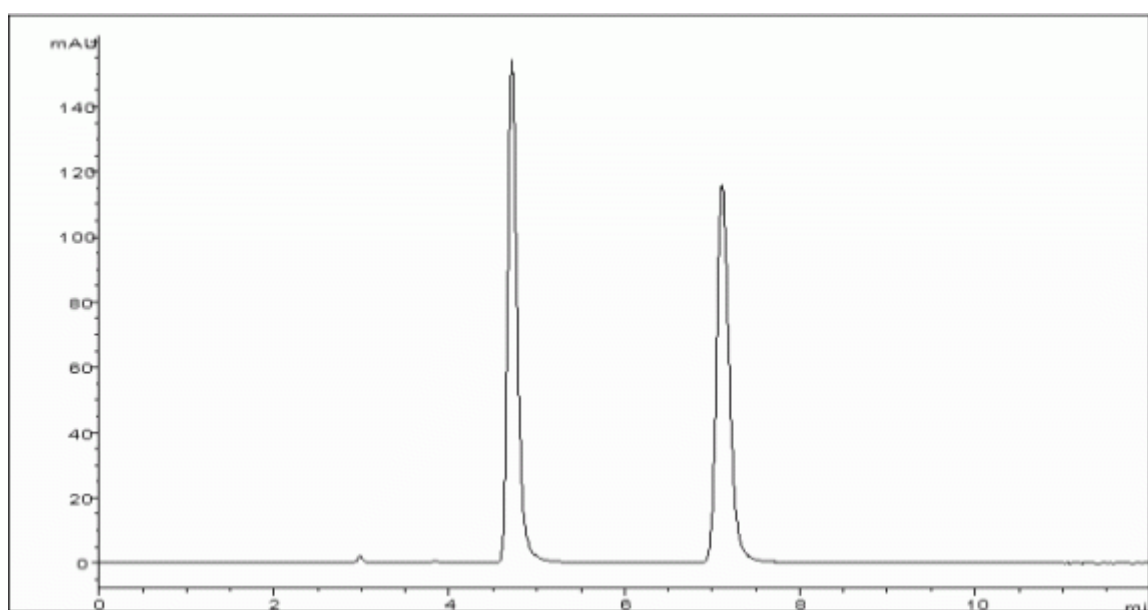

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethyl acetate / ethanol / ethanalamine = 50 / 50 / 2 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.006 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.7	0.57	1 - 2	10.96
2	7.1	1.37		