


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Phenylsuccinic acid

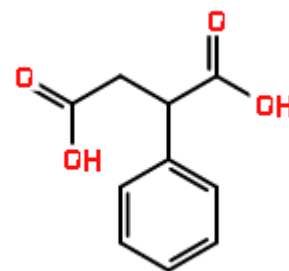

CHROMATOGRAM

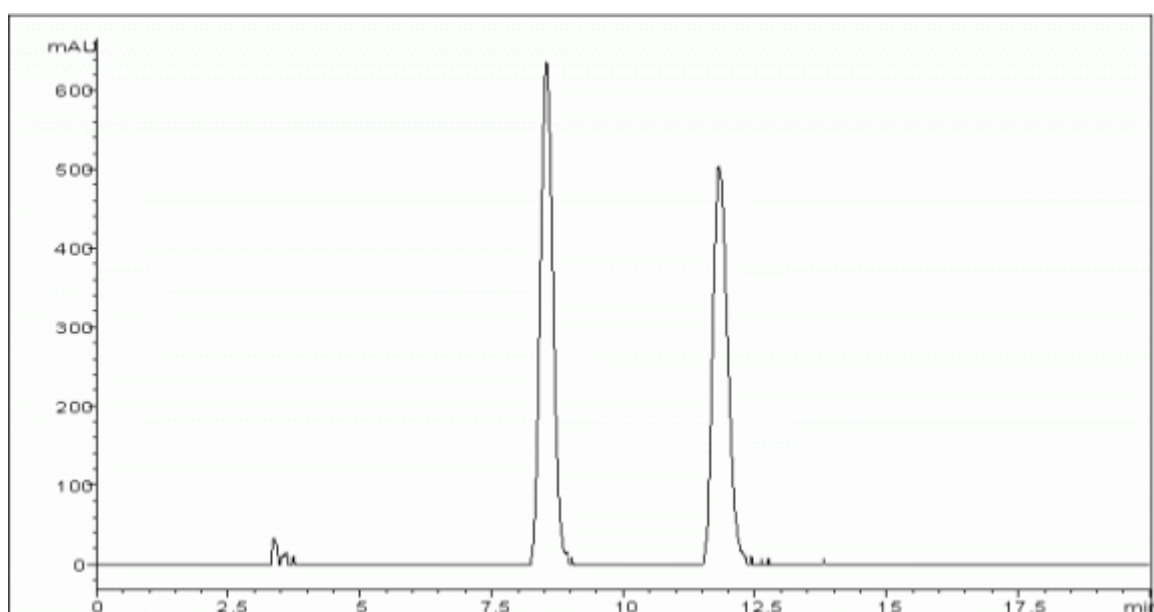

ANALYTICAL REPORT

| | | | |
|-----------------------------|--|----------------------|----|
| Column Trade Name | CHIRALPAK® IA | | |
| Column Size (Particle Size) | 4.6 x 250 mm (5 µm) | | |
| Mobile Phase Composition | dichloromethane / methanol / trifluoroacetic acid = 98 / 2 / 0.1 | | |
| Chromatographic Mode | Normal Phase | | |
| Flow rate (mL/min) | 1 | Temperature (°C) | 25 |
| Injection Volume (µL) | 0.02 (mg) | Sample Concentration | |
| Detection | ELSD | | |
| Elution Order | | | |

| No. | Rt | K' | Alpha Value | Resolution |
|-----|------|------|-------------|------------|
| 1 | 8.5 | 1.80 | 1 - 2 | 6.71 |
| 2 | 11.8 | 2.88 | | |