


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Sertraline (enriched)

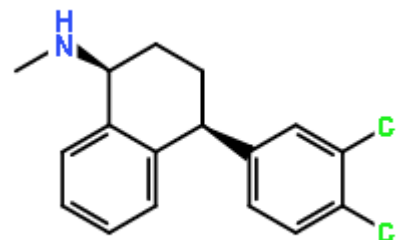

CHROMATOGRAM

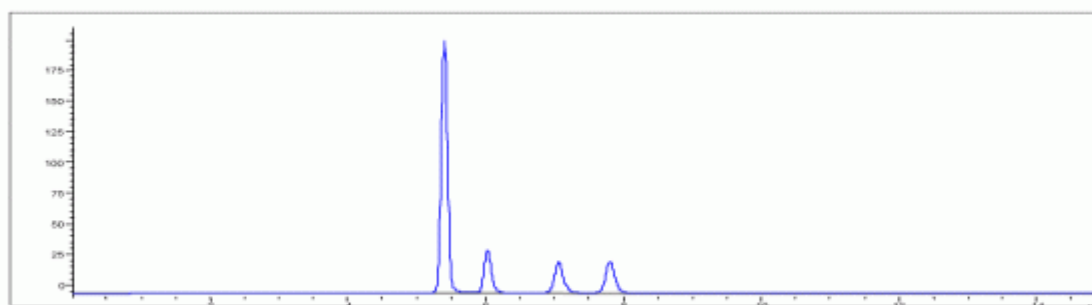

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / methanol / diethylamine = 98 / 1 / 1 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.4	0.81	1 - 2	3.60
2	6.0	1.01	2 - 3	5.00
3	7.1	1.35	3 - 4	3.20
4	7.8	1.60		