


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Tirofiban (racemic)

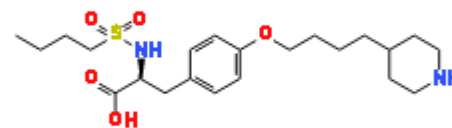

CHROMATOGRAM

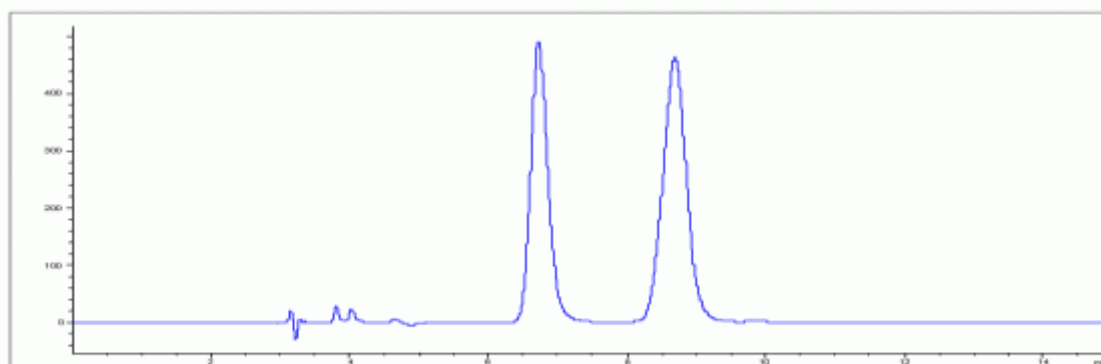

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / methanol / 2-propanol / trifluoroacetic acid / ethylenediamine = 30 / 30 / 40 / 0.3 / 0.3		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.7	1.24	1 - 2	3.51
2	8.7	1.91		