


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Tropic acid

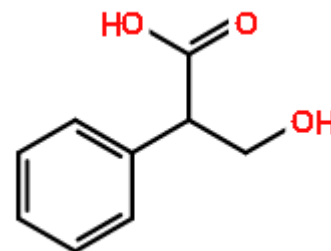

## CHROMATOGRAM

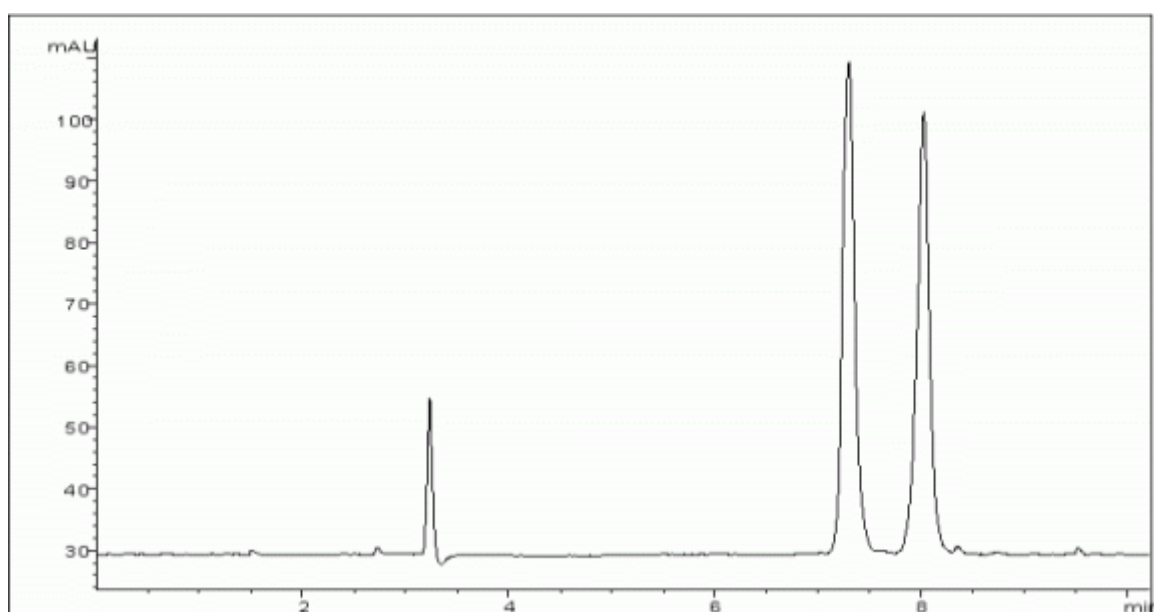

## ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IC 4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethyl acetate / trifluoroacetic acid = 70 / 30 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.3	1.29	1 - 2 1.17	3.51
2	8.0	1.51		