


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

3-oxo-1-Indancarboxylic acid

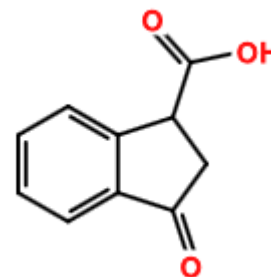

CHROMATOGRAM

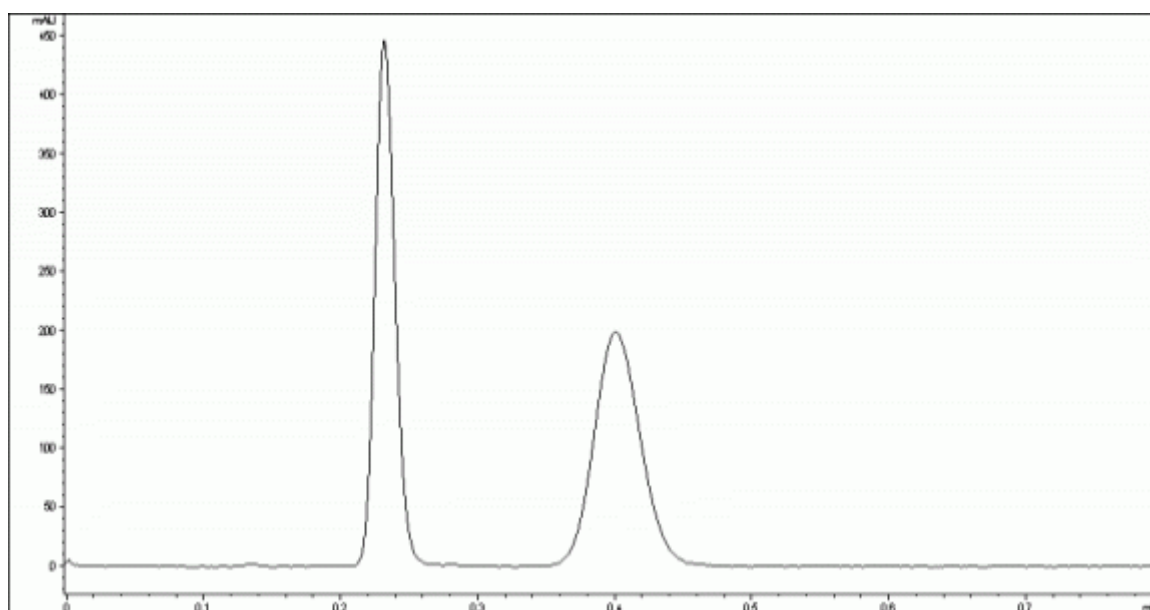

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AY-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / trifluoroacetic acid = 70 / 30 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	5	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.2	0.91	1 - 2	2.52
2	0.4	2.31		
				3.83