


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Chlorthalidone

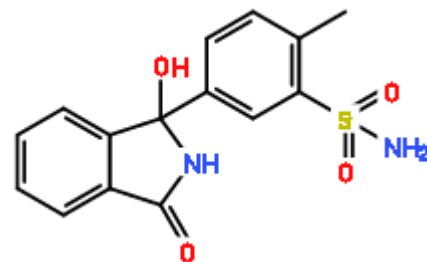

CHROMATOGRAM

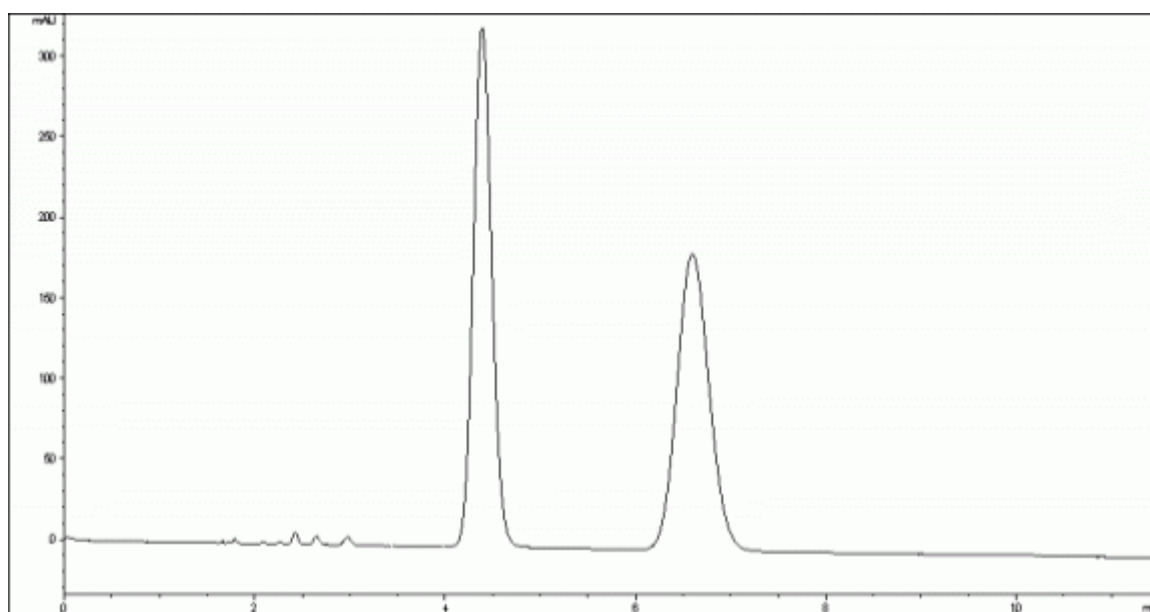

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL [®] OZ-3		
Column Size (Particle Size)	4.6 x 150 mm (3 μm)		
Mobile Phase Composition	n-hexane / ethanol = 60 / 40		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (μL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.4	1.60	1 - 2	4.51
2	6.6	2.90		