


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Mecoprop

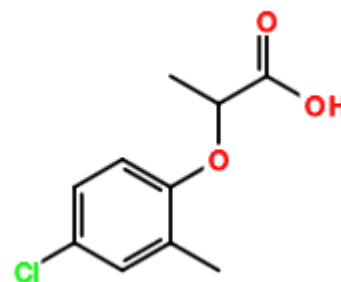

CHROMATOGRAM

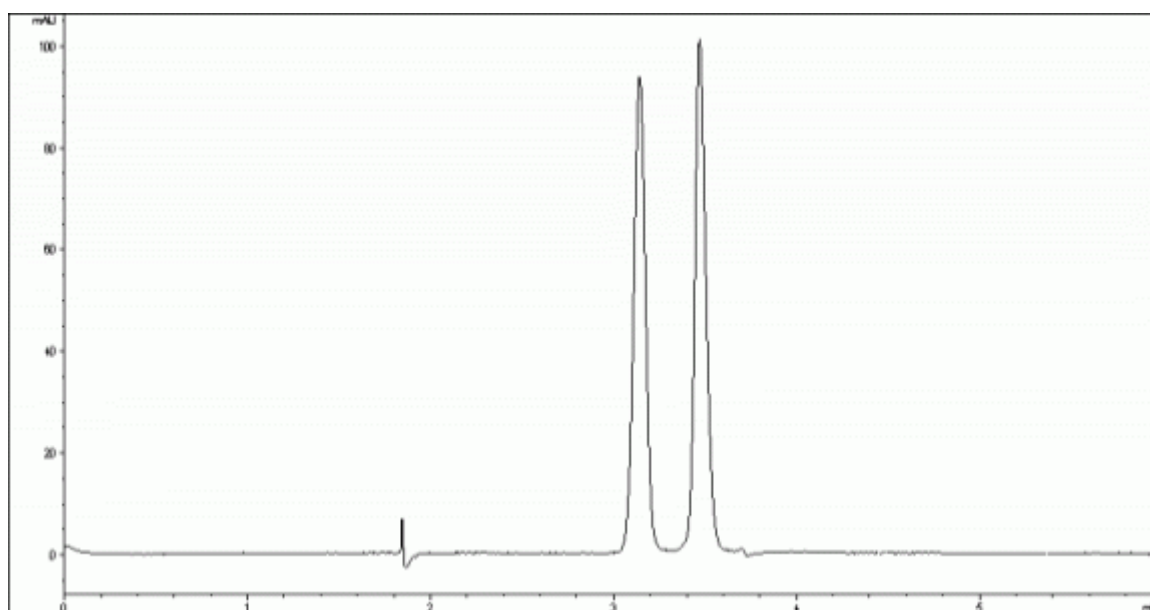

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® AY-3 4.6 x 150 mm (3 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / trifluoroacetic acid = 95 / 5 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.1	0.78	1 - 2	1.24
2	3.5	0.96		