


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Methaqualone

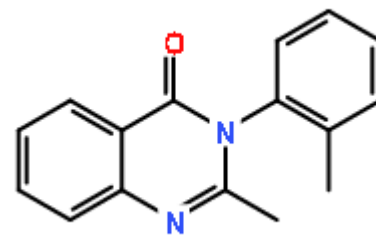

CHROMATOGRAM

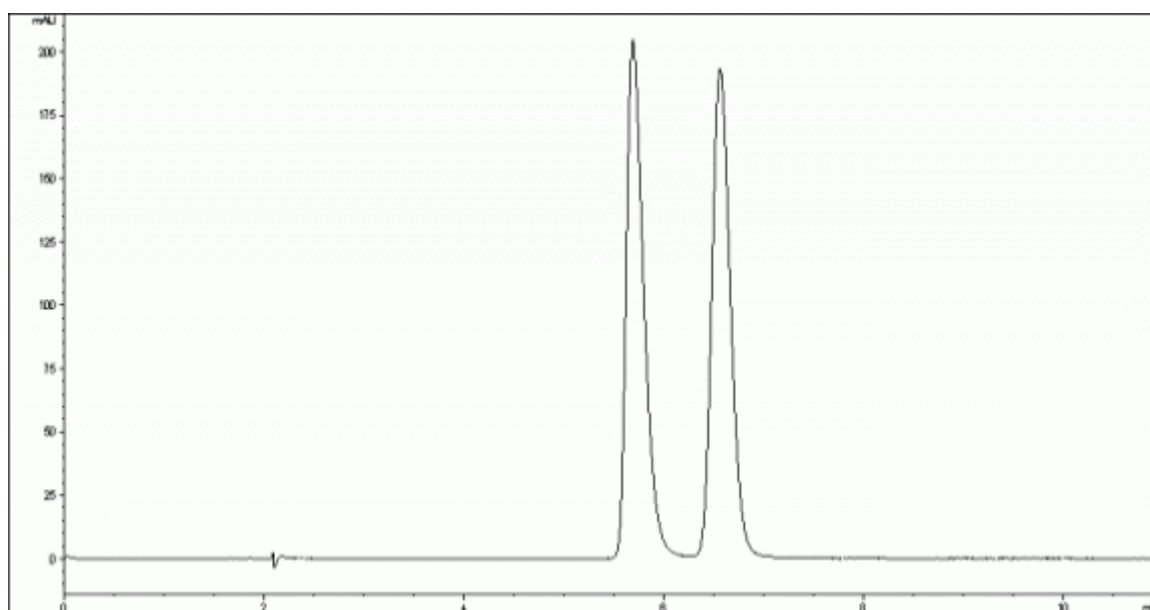

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL [®] OZ-3		
Column Size (Particle Size)	4.6 x 150 mm (3 μ m)		
Mobile Phase Composition	n-hexane / 2-propanol = 80 / 20		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature ($^{\circ}$ C)	25
Injection Volume (μ L)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.7	2.37	1 - 2	2.68
2	6.6	2.88		