


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Econazole

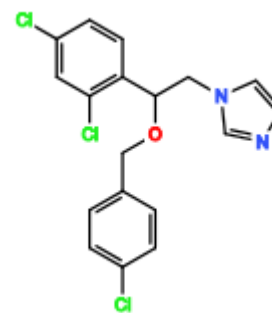

CHROMATOGRAM

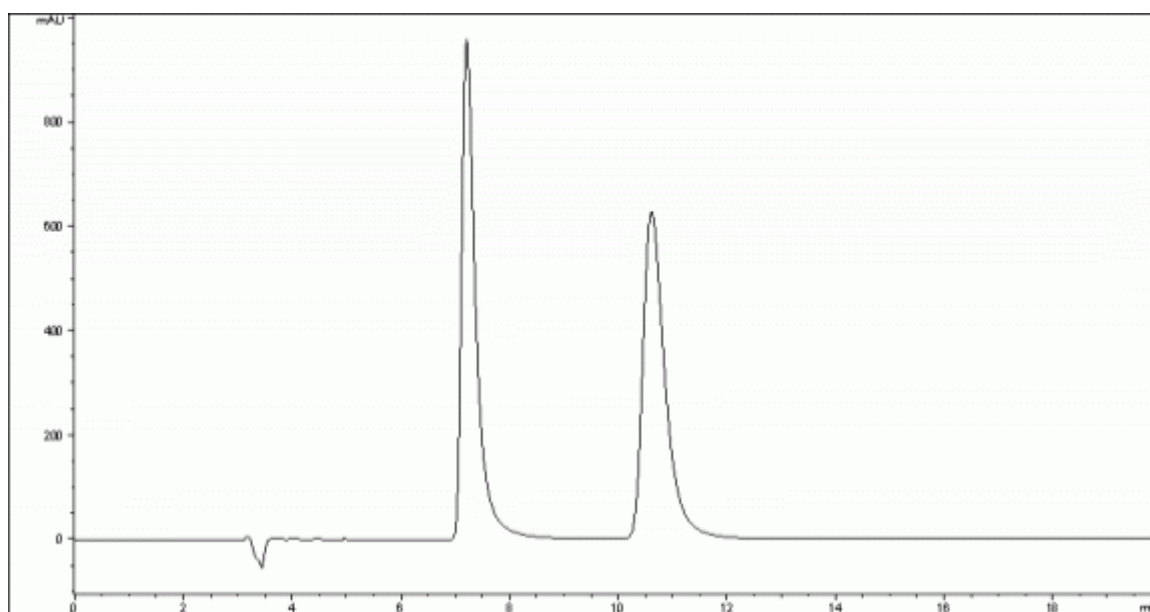

ANALYTICAL REPORT

| | | | |
|-----------------------------|--|----------------------|----|
| Column Trade Name | CHIRALCEL® OZ-H | | |
| Column Size (Particle Size) | 4.6 x 250 mm (5 µm) | | |
| Mobile Phase Composition | n-hexane / 2-propanol / diethylamine = 50 / 50 / 0.1 | | |
| Chromatographic Mode | Normal Phase | | |
| Flow rate (mL/min) | 1 | Temperature (°C) | 25 |
| Injection Volume (µL) | 0.01 (mg) | Sample Concentration | |
| Detection | UV-VIS 230 nm | | |
| Elution Order | | | |

| No. | Rt | K' | Alpha Value | Resolution |
|-----|------|------|-------------|------------|
| 1 | 7.2 | 1.40 | 1 - 2 | 6.00 |
| 2 | 10.6 | 2.53 | | |