


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Chlorpheniramine

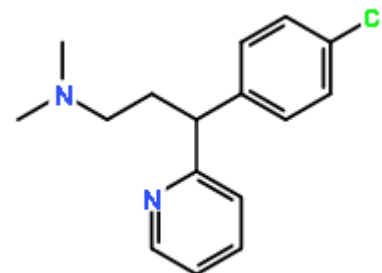

## CHROMATOGRAM

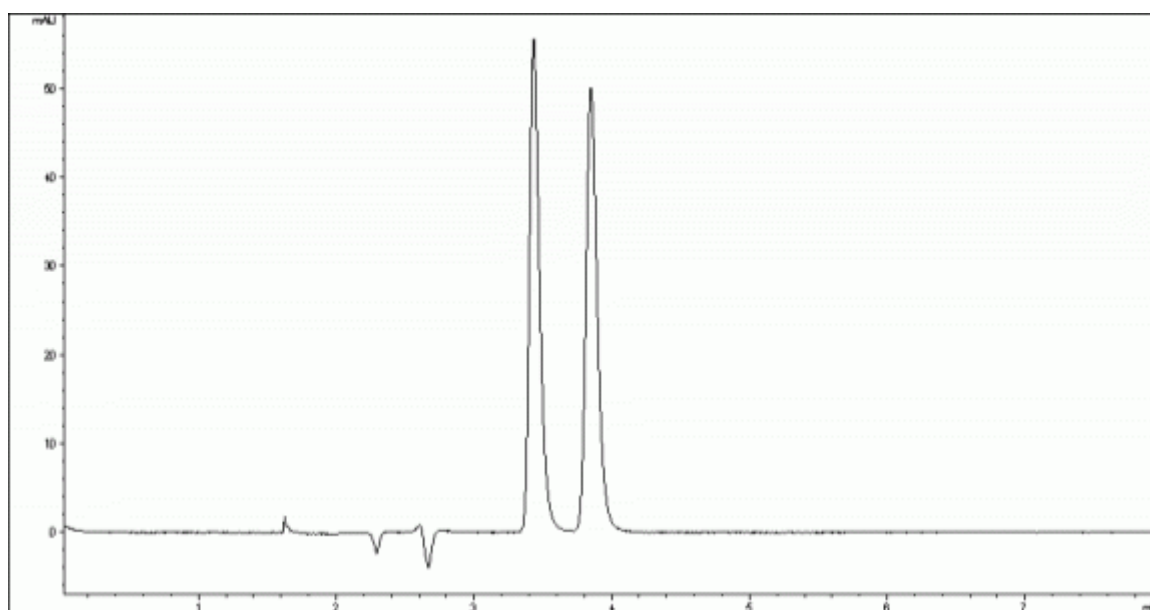

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / diethylamine = 90 / 10 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.001 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.4	1.12	1 - 2	3.03
2	3.9	1.38		