


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Dorzolamide (racemic)

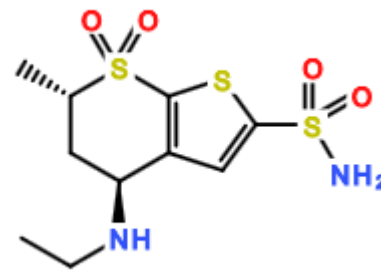

CHROMATOGRAM

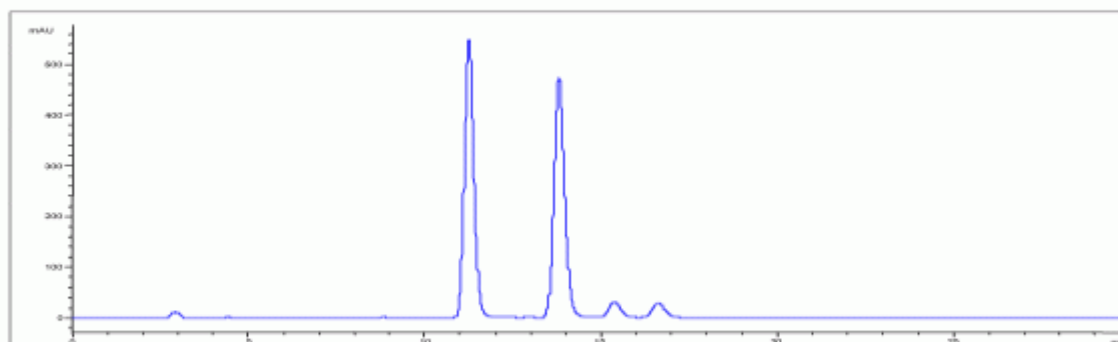

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 80 / 20 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	11.3	2.77	1 - 2	4.62
2	13.8	3.61	2 - 3	2.64
3	15.4	4.14	3 - 4	1.91
4	16.6	4.54		