


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Lorazepam

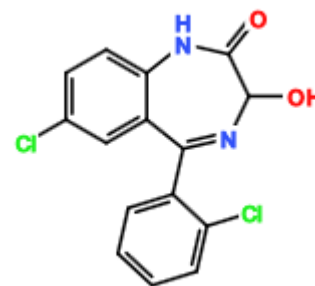

CHROMATOGRAM

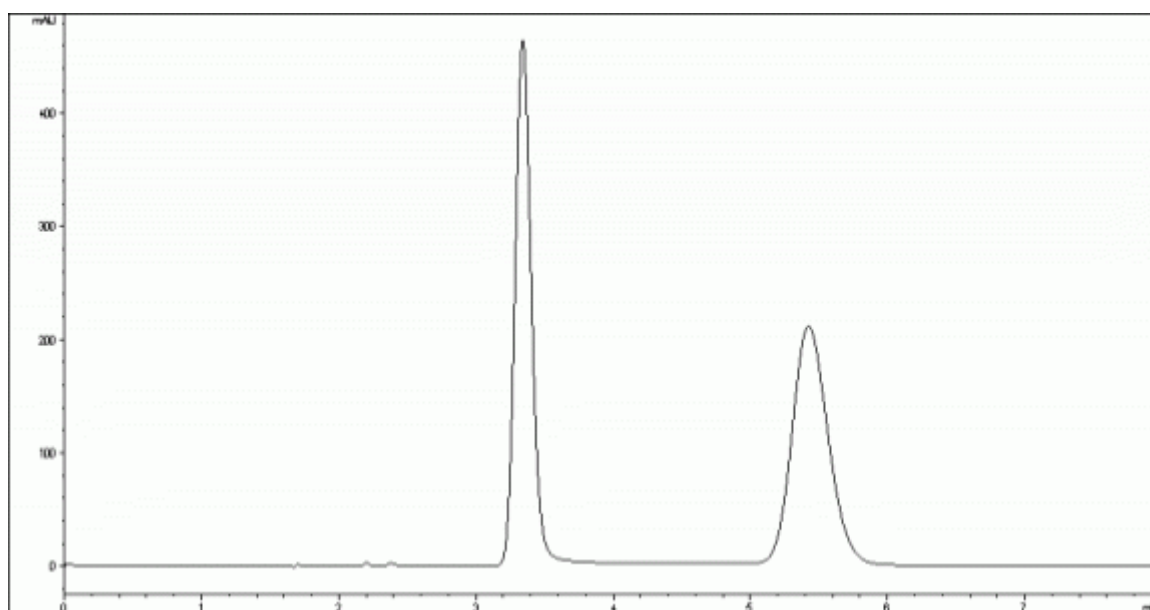

ANALYTICAL REPORT

| | | | |
|-----------------------------|------------------------------|----------------------|----|
| Column Trade Name | CHIRALPAK® AS-3 | | |
| Column Size (Particle Size) | 4.6 x 150 mm (3 µm) | | |
| Mobile Phase Composition | n-hexane / ethanol = 50 / 50 | | |
| Chromatographic Mode | Normal Phase | | |
| Flow rate (mL/min) | 1 | Temperature (°C) | 25 |
| Injection Volume (µL) | 0.005 (mg) | Sample Concentration | |
| Detection | UV-VIS 254 nm | | |
| Elution Order | | | |

| No. | Rt | K' | Alpha Value | Resolution |
|-----|-----|------|-------------|------------|
| 1 | 3.3 | 0.99 | 1 - 2 | 6.02 |
| 2 | 5.4 | 2.23 | | |