


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Abscisic acid

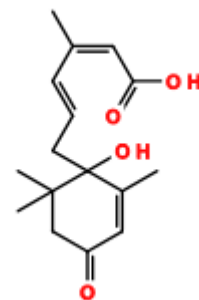

## CHROMATOGRAM

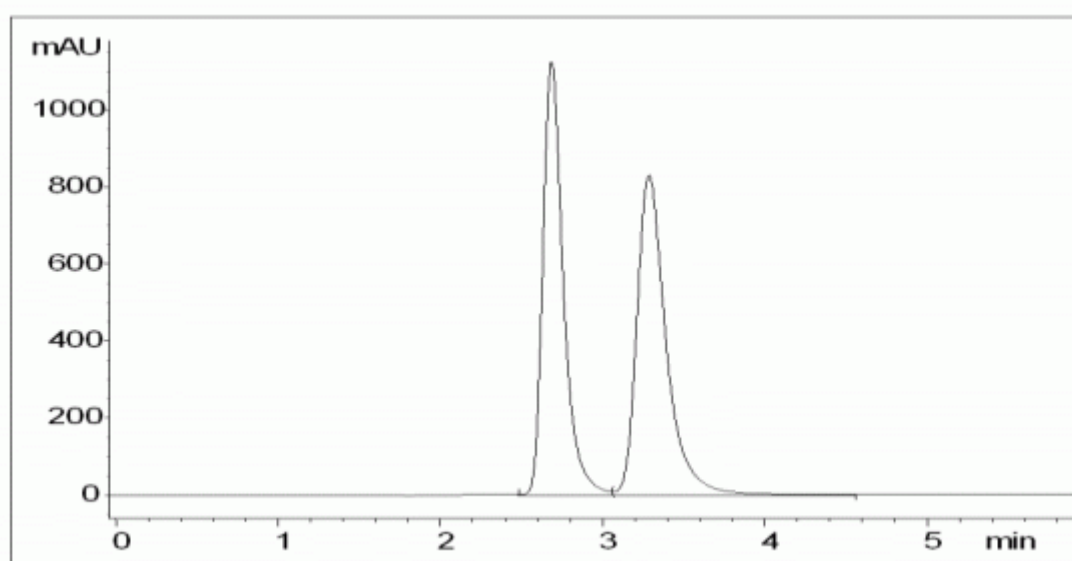

## ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OD-H		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	n-heptane / 2-propanol / trifluoroacetic acid = 80 / 20 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.7	0.48	1 - 2	2.90
2	3.3	0.83		