


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Citalopram

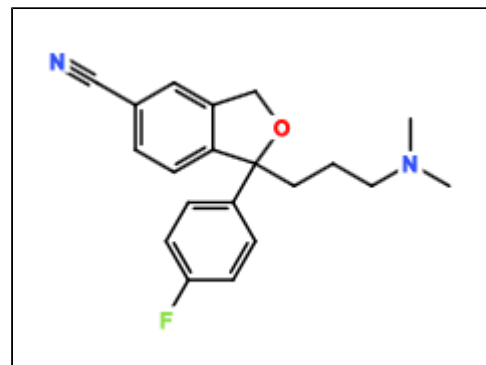

CHROMATOGRAM

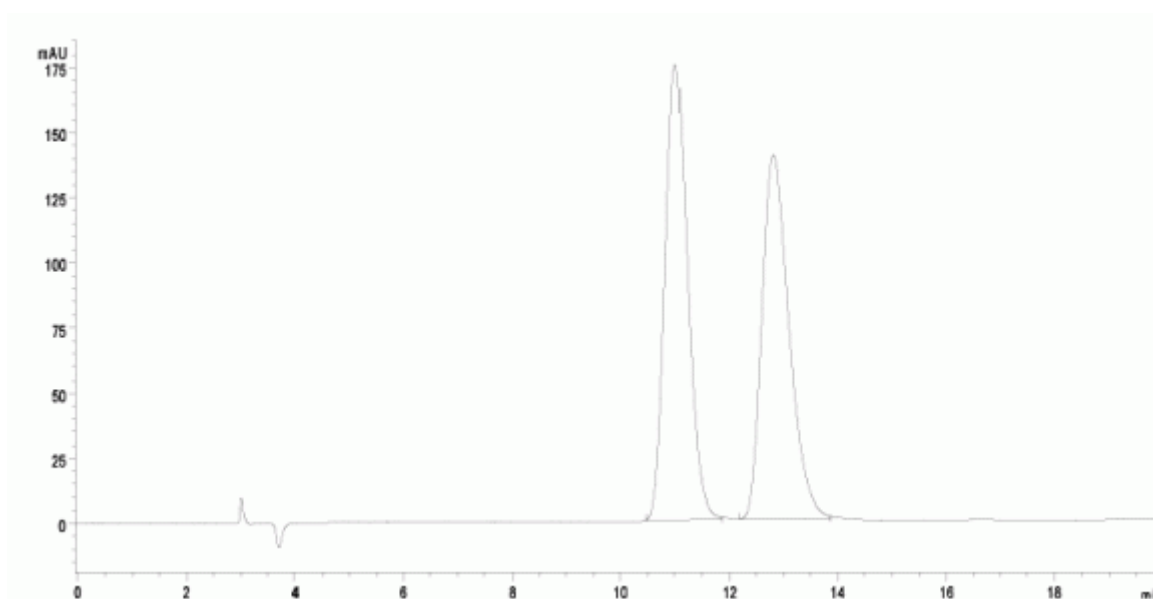

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL [®] OD-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-heptane / 2-propanol / diethylamine = 94 / 6 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 250 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	11.0	2.66	1 - 2	2.20
2	12.8	3.27		