


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

DL-Leu-DL-Val

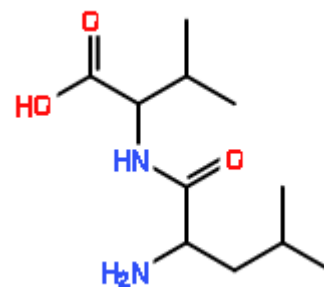

CHROMATOGRAM

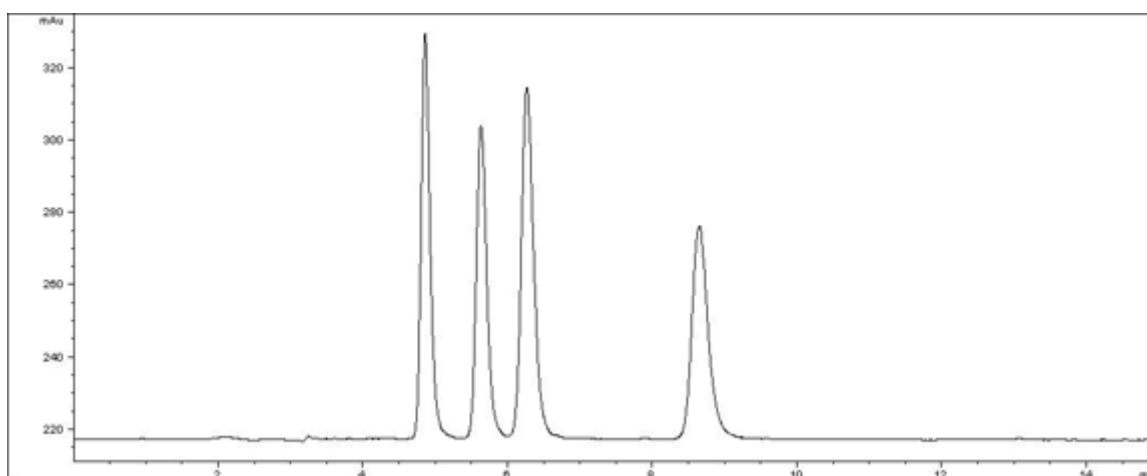

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ZWIX(+) 3 x 150 mm (3 µm)		
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / tetrahydrofuran / water = 49 / 49 / 2		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.9	1.72	1 - 2	1.25
2	5.6	2.15	2 - 3	1.17
3	6.3	2.50	3 - 4	1.53
4	8.7	3.84		6.81