


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Tolterodine (racemic)

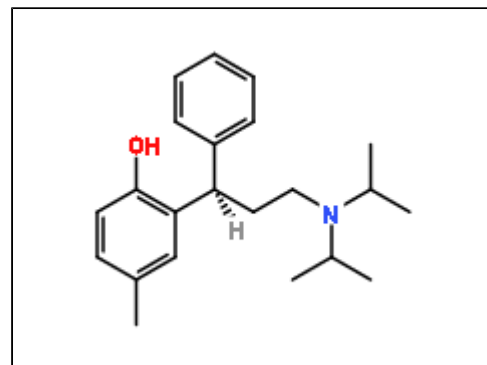

CHROMATOGRAM

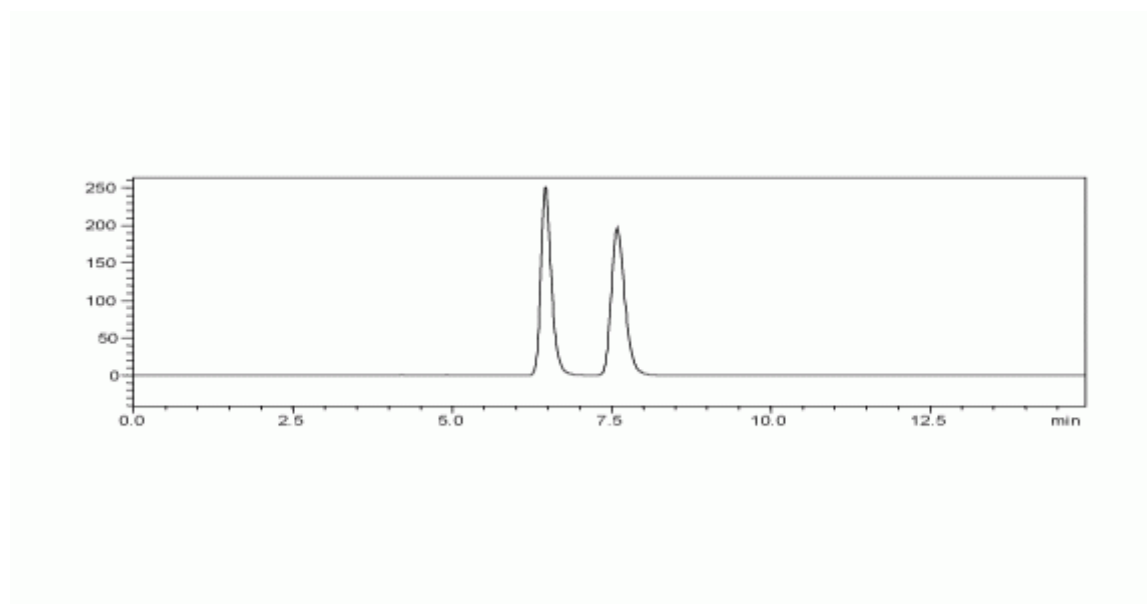

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AS-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / 2-butanol / diethylamine = 98 / 1 / 1 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	0.8	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 285 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.4	0.72	1 - 2	3.00
2	7.6	1.03		