


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

4-Phenyl-2-oxazolidinone

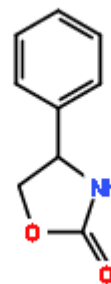

## CHROMATOGRAM

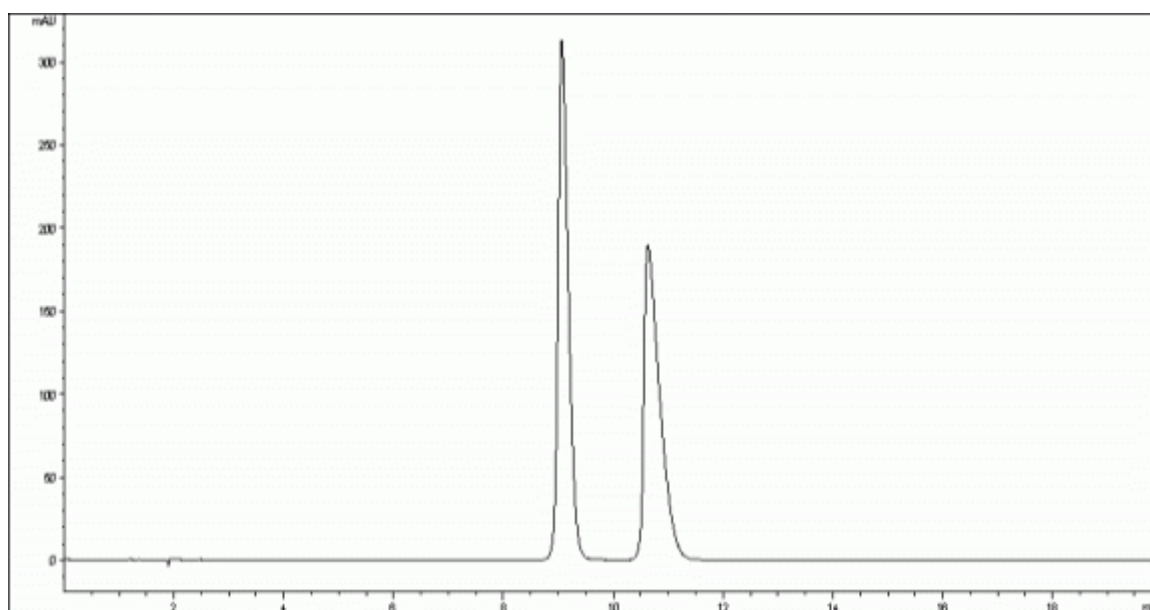

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	water / acetonitrile = 80 / 20		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.1	4.19	1 - 2	3.72
2	10.6	5.08		