


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Ornidazole

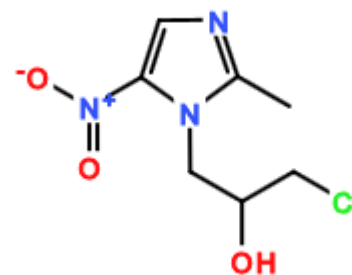

CHROMATOGRAM

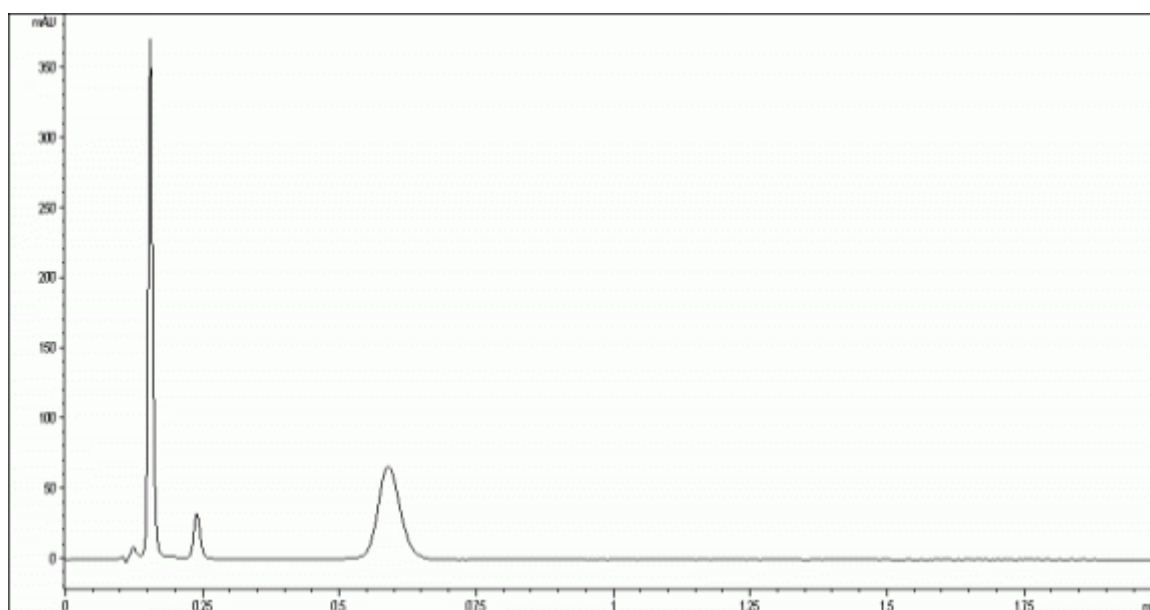

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	20mM NH ₄ CO ₃ aq. pH9.0 / acetonitrile = 30 / 70		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	5	Temperature (°C)	40
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.2	0.30	1 - 2	13.08
2	0.6	3.96		9.75