


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Bicalutamide

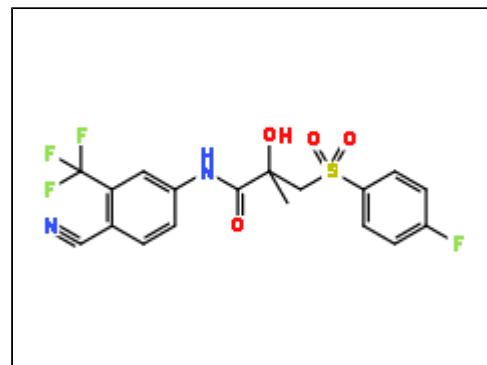

CHROMATOGRAM

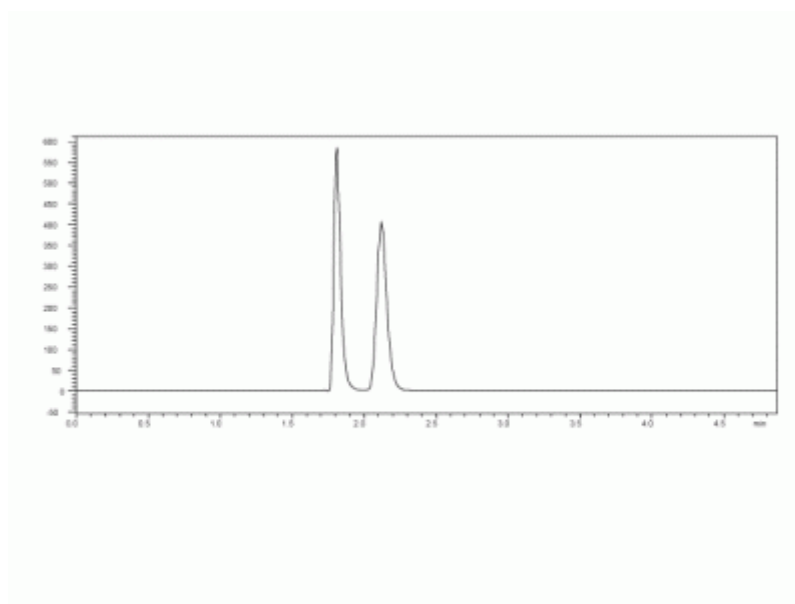

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-3R		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	5mM NH ₄ OAc aq. / methanol / acetonitrile = 5 / 60 / 35		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	1.8	0.01	1 - 2	17.00
2	2.1	0.19		2.53