


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Flavanone

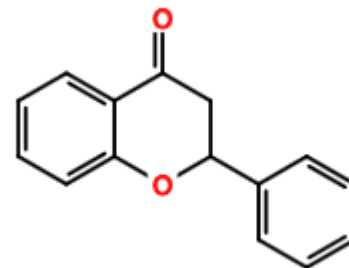

## CHROMATOGRAM

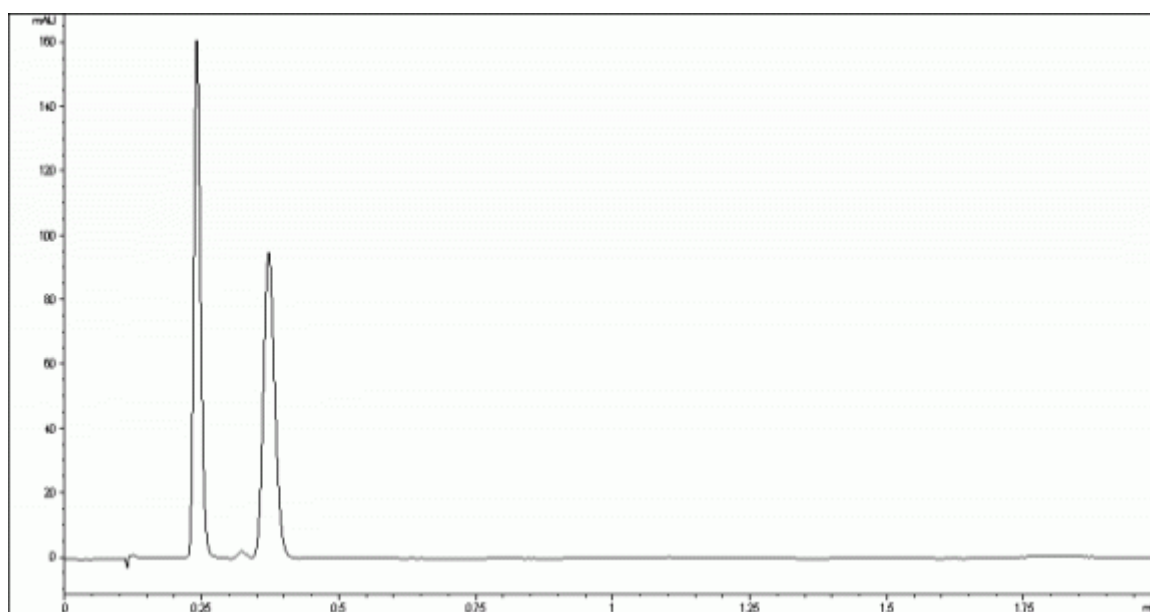

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-3R		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	water / acetonitrile = 20 / 80		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	5	Temperature (°C)	40
Injection Volume (µL)	0.001 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.2	0.96	1 - 2	4.38
2	0.4	2.02		