


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Tolperisone

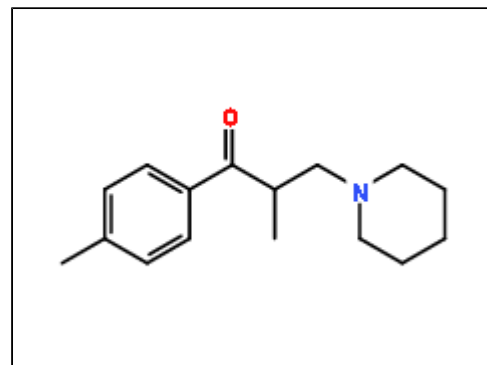

CHROMATOGRAM

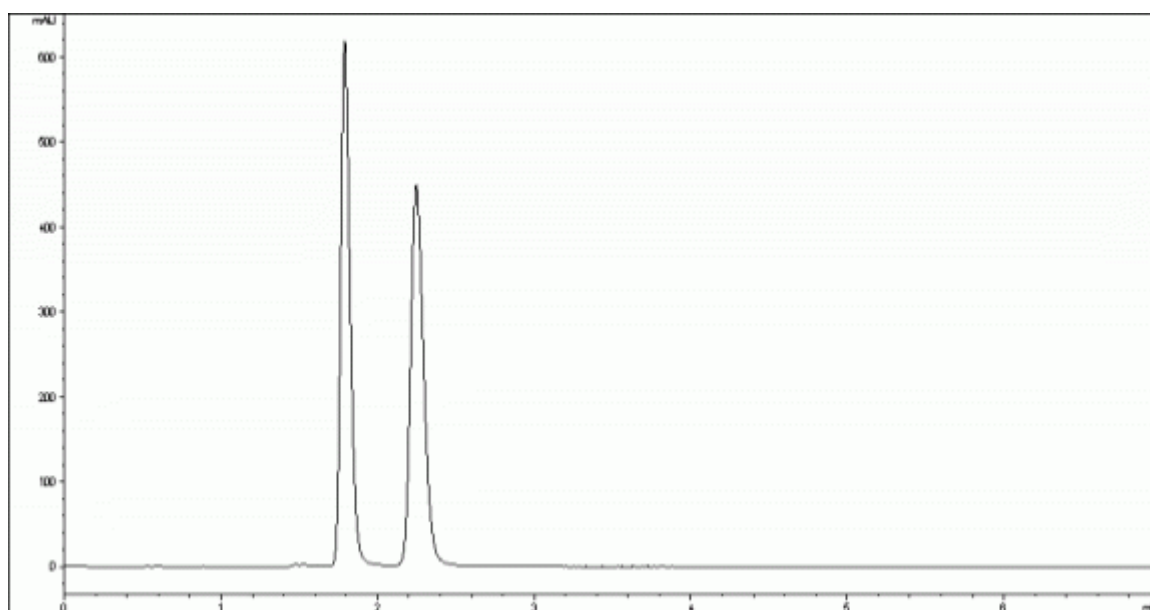

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-3R		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	20mM NH ₄ CO ₃ aq. pH9.0 / acetonitrile = 40 / 60		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	1.8	1.91	1 - 2	3.55
2	2.3	2.66		