


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

alpha-Methyl-alpha-phenylsuccinimide

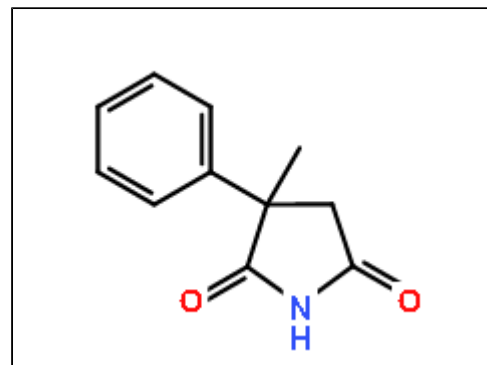

CHROMATOGRAM

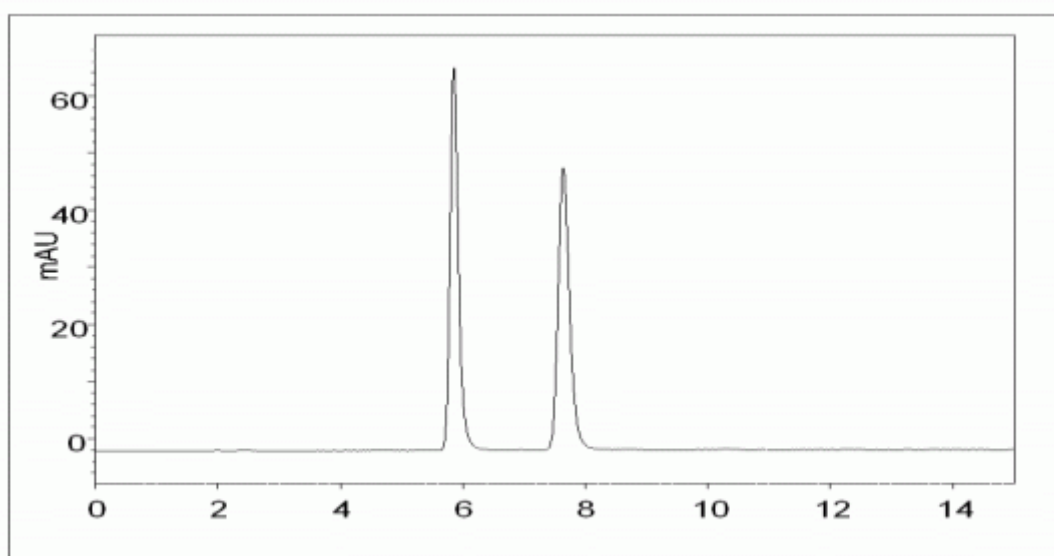

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AS-RH		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	water / acetonitrile = 50 / 50		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.9	0.63	1 - 2	9.10
2	7.6	1.12		