


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Methyl-(phenylsulfinyl)-acetate

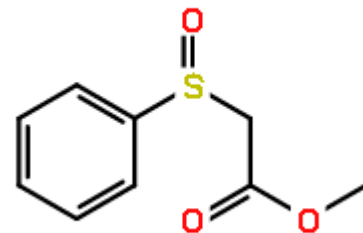

## CHROMATOGRAM

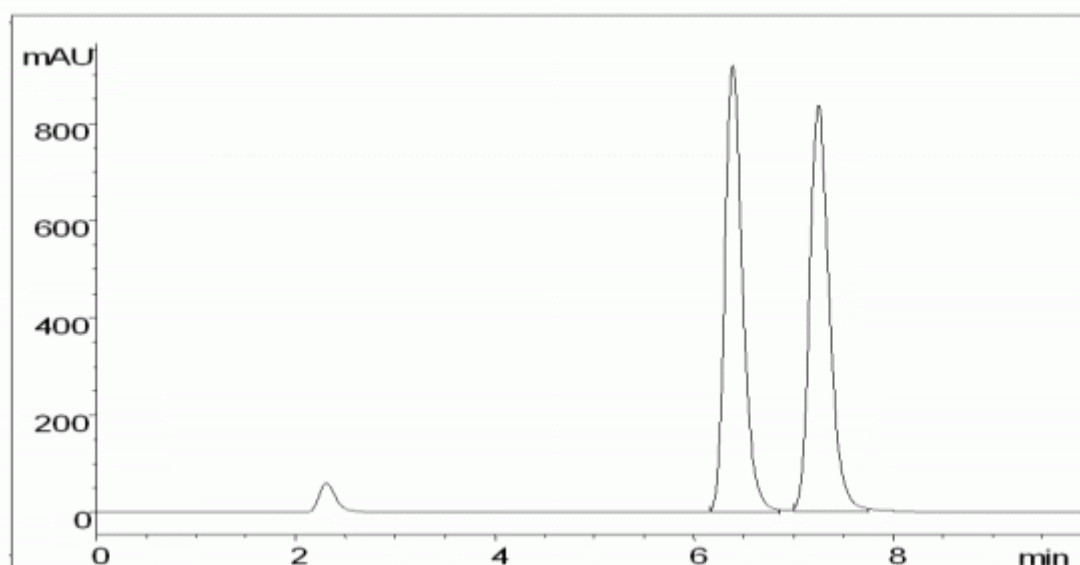

## ANALYTICAL REPORT

|                             |                                |                      |    |
|-----------------------------|--------------------------------|----------------------|----|
| Column Trade Name           | CHIRALPAK® AS-RH               |                      |    |
| Column Size (Particle Size) | 4.6 x 150 mm (5 µm)            |                      |    |
| Mobile Phase Composition    | water / acetonitrile = 60 / 40 |                      |    |
| Chromatographic Mode        | Reverse Phase                  |                      |    |
| Flow rate (mL/min)          | 0.5                            | Temperature (°C)     | 25 |
| Injection Volume (µL)       | 0.01 (mg)                      | Sample Concentration |    |
| Detection                   | UV-VIS 230 nm                  |                      |    |
| Elution Order               |                                |                      |    |

| No. | Rt  | K'   | Alpha Value | Resolution |
|-----|-----|------|-------------|------------|
| 1   | 6.4 | 0.60 | 1 - 2       | 2.50       |
| 2   | 7.3 | 0.80 |             |            |