


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Troger's base

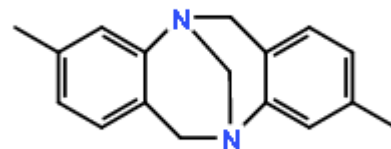

CHROMATOGRAM

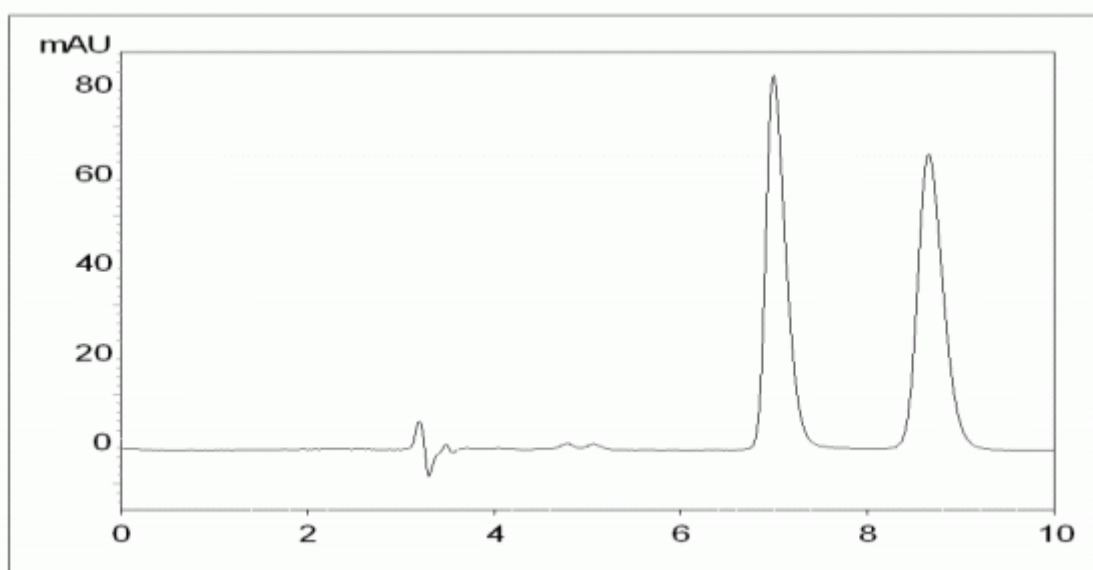

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AS-RH		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	H3PO4 aq. pH2.0 / acetonitrile = 65 / 35		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.0	2.88	1 - 2	3.71
2	8.6	3.80		