


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Ketoprofen

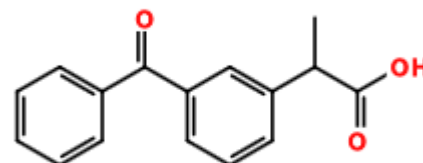

CHROMATOGRAM

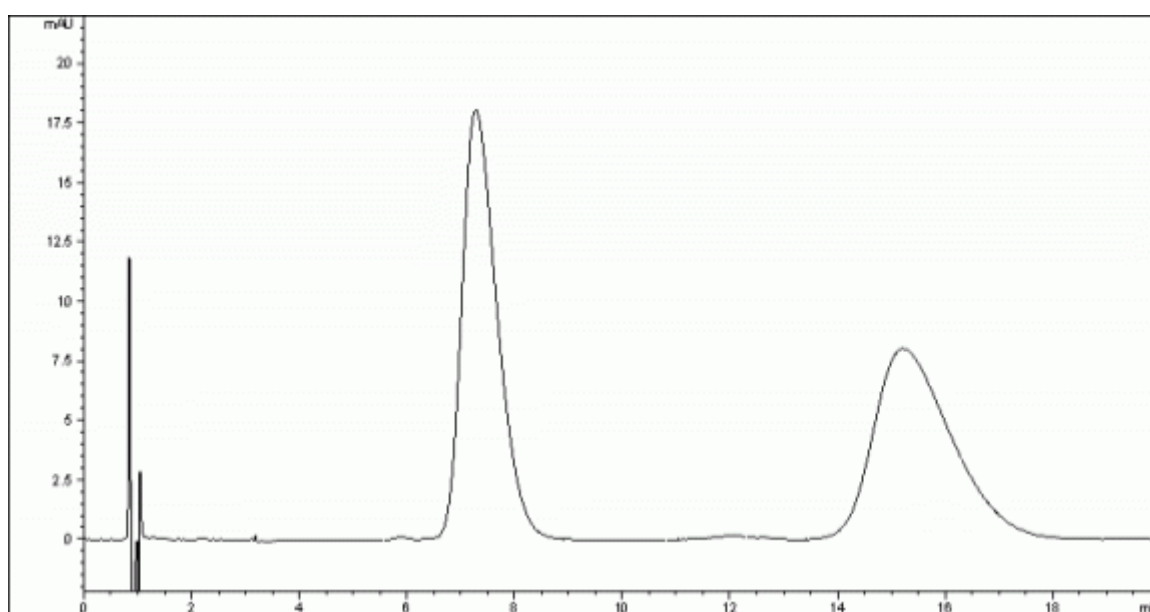

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AGP		
Column Size (Particle Size)	4 x 100 mm (5 µm)		
Mobile Phase Composition	10mM NH ₄ OAc aq. (+5mM N,N-dimethyloctyl amine) pH4.5 / 2-propanol = 98 / 2		
Chromatographic Mode	Reverse Phase		

Flow rate (mL/min)	0.9	Temperature (°C)	20
Injection Volume (µL)	0.0005 (mg)	Sample Concentration	
Detection	UV-VIS 225 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.3	5.88	1 - 2	4.14
2	15.2	13.37		