


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Ketorolac

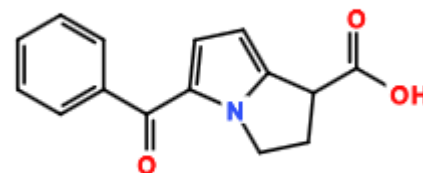

CHROMATOGRAM

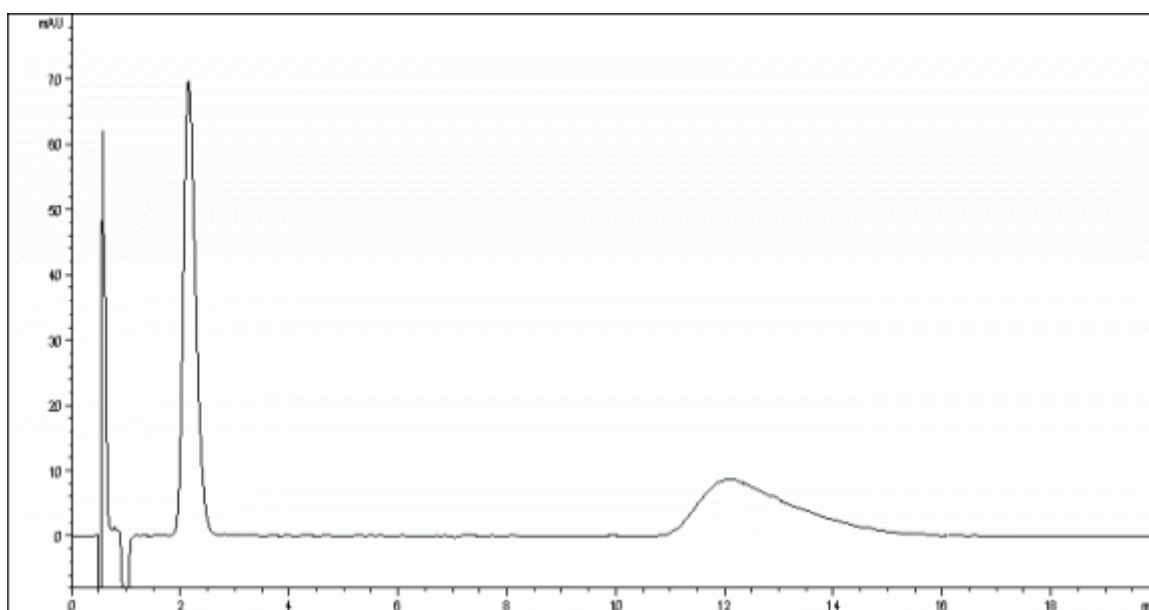

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® HSA		
Column Size (Particle Size)	4 x 50 mm (5 µm)		
Mobile Phase Composition	10mM NH ₄ OAc aq. (+2mM octanoic acid aq.) pH7.0 / 2-propanol = 85 / 15		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.9	Temperature (°C)	30
Injection Volume (µL)	0.0005 (mg)	Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.2	2.74	1 - 2	5.37
2	12.1	20.10		