


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Gly-DL-Leu-DL-Ala

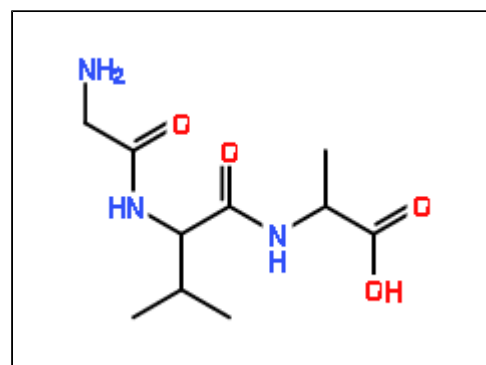

CHROMATOGRAM

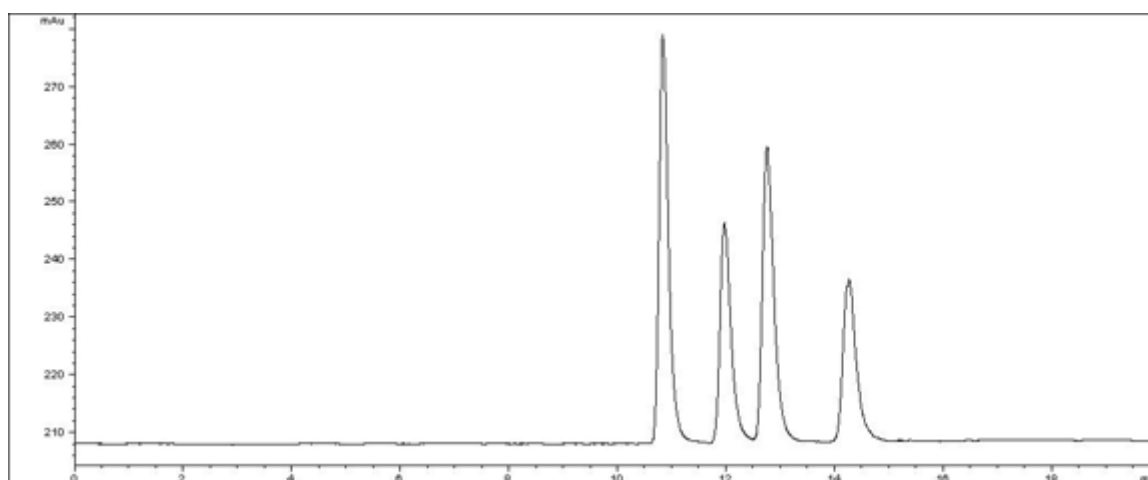

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	4 x 250 mm (3 µm)		
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / water = 98 / 2		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	10
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	10.8	1.15	1 - 2	3.48
2	12.0	1.38	2 - 3	2.17
3	12.8	1.54	3 - 4	3.79
4	14.3	1.84		