


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Sotalol

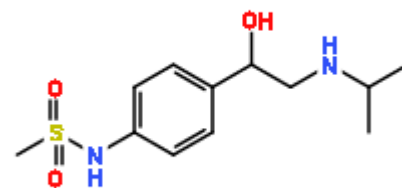

CHROMATOGRAM

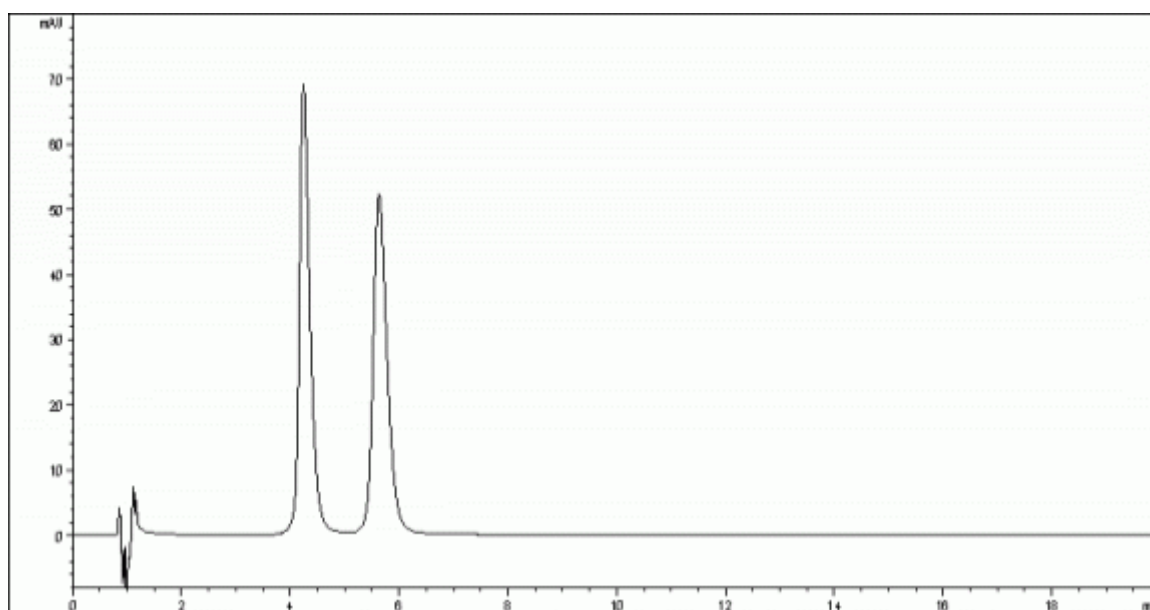

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® CBH		
Column Size (Particle Size)	4 x 100 mm (5 µm)		
Mobile Phase Composition	10mM NH ₄ OAc aq. pH7.0 / methanol = 95 / 5		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.9	Temperature (°C)	20
Injection Volume (µL)	0.0005 (mg)	Sample Concentration	
Detection	UV-VIS 225 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.2	2.98	1 - 2	3.54
2	5.6	4.29		