


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Gly-DL-Val

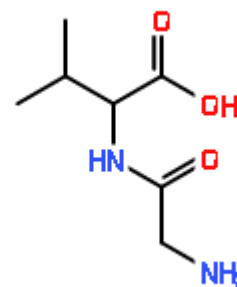

CHROMATOGRAM

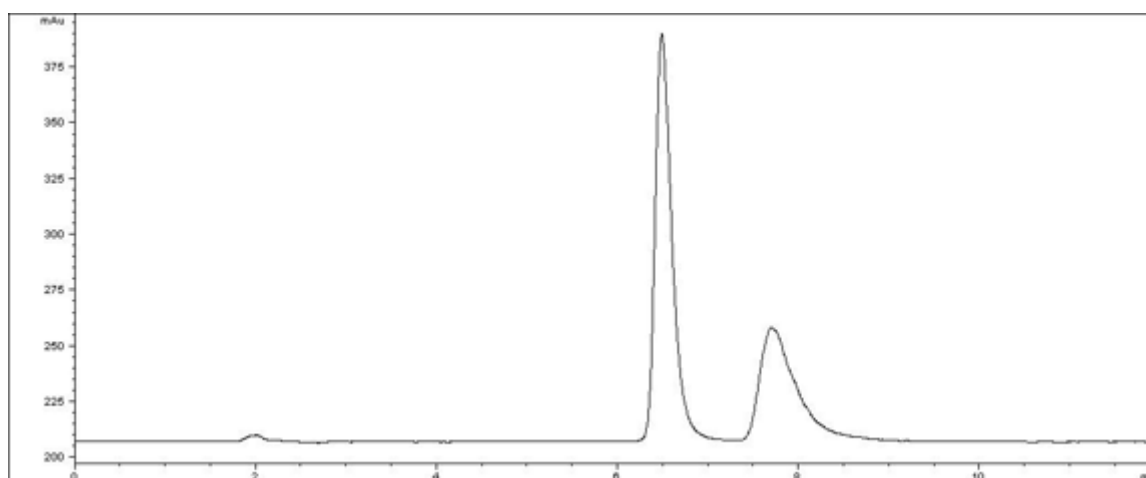

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	3 x 150 mm (3 µm)		
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / acetonitrile / water = 49 / 49 / 2		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.5	2.65	1 - 2	2.35
2	7.7	3.33		