

## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Tetramisole



## CHROMATOGRAM

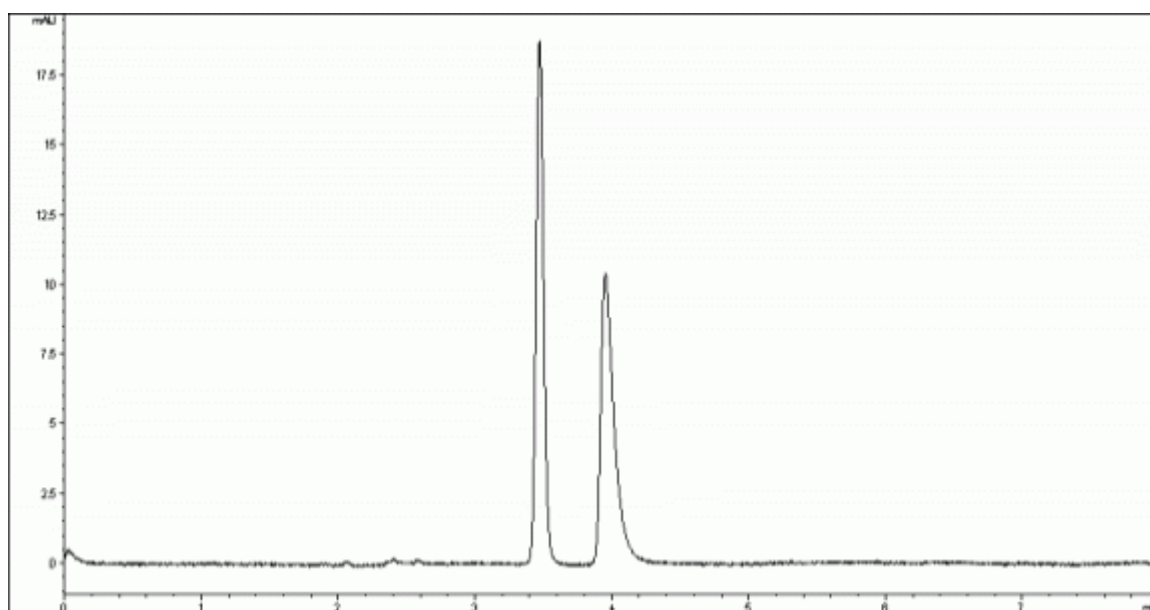

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	acetonitrile / ethanolamine = 100 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.001 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.5	1.04	1 - 2	3.61
2	4.0	1.32		