


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Thalidomide

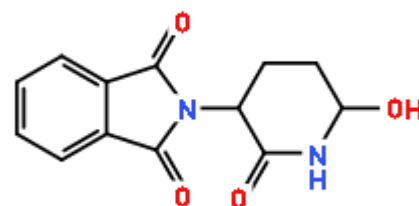

CHROMATOGRAM

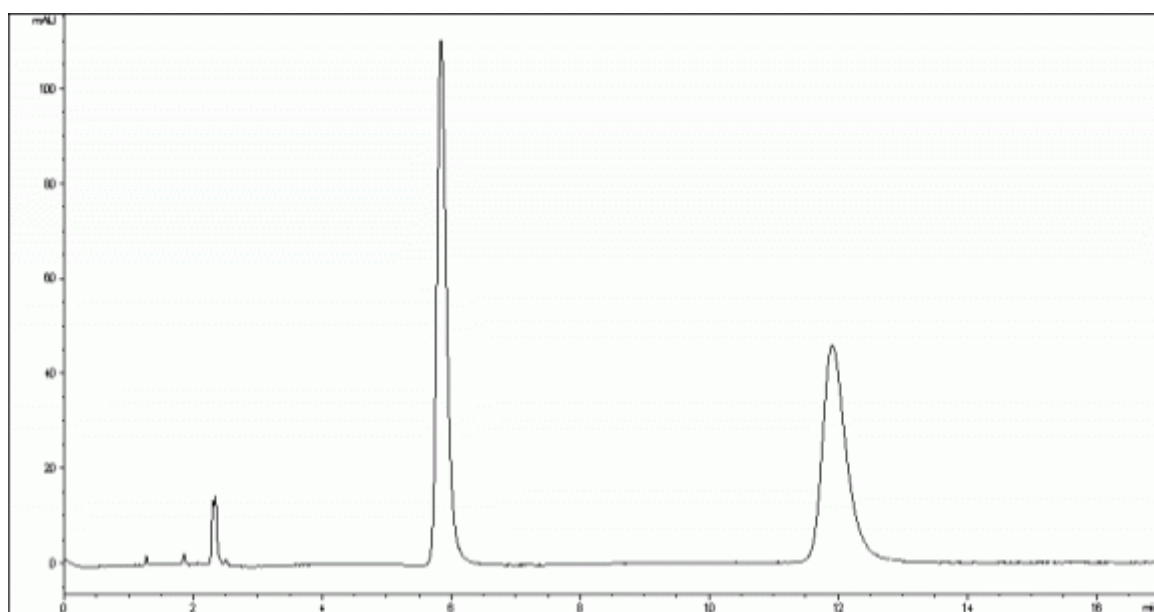

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	methanol = 100		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.001 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.8	2.43	1 - 2	13.32
2	11.9	5.99		