


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Kavain

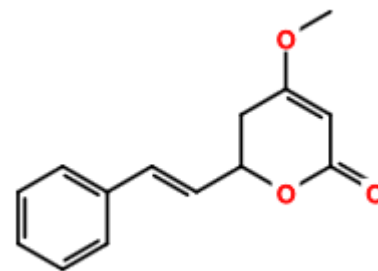

CHROMATOGRAM

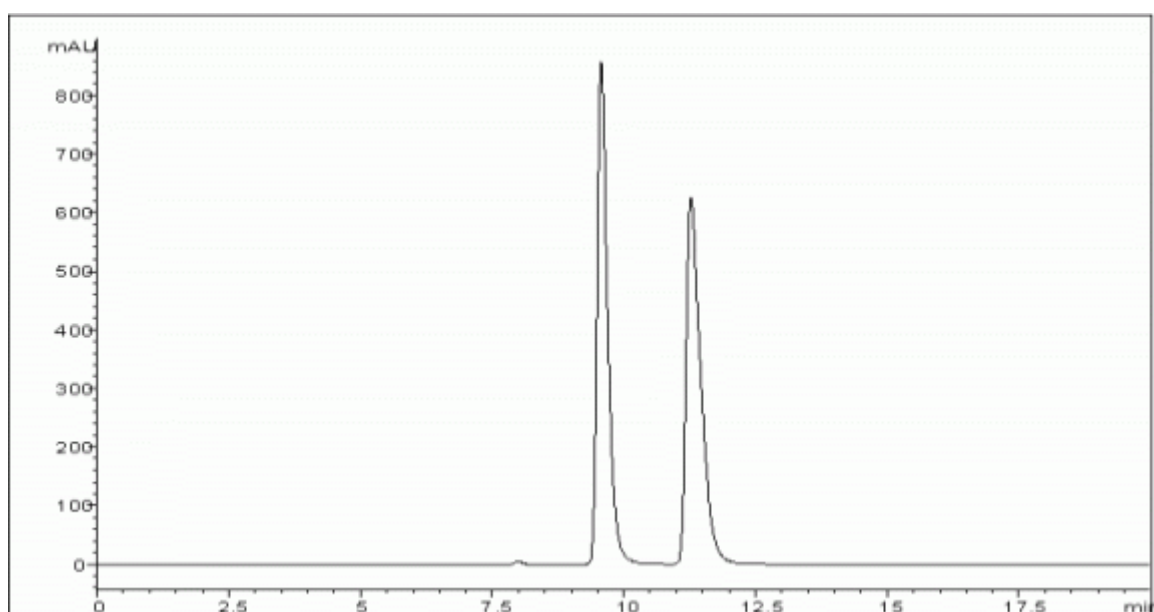

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	methanol = 100		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.7	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.6	1.29	1 - 2	4.06
2	11.3	1.70		