


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Pramipexole (enriched)

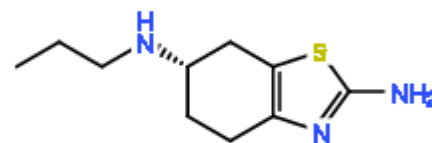

## CHROMATOGRAM

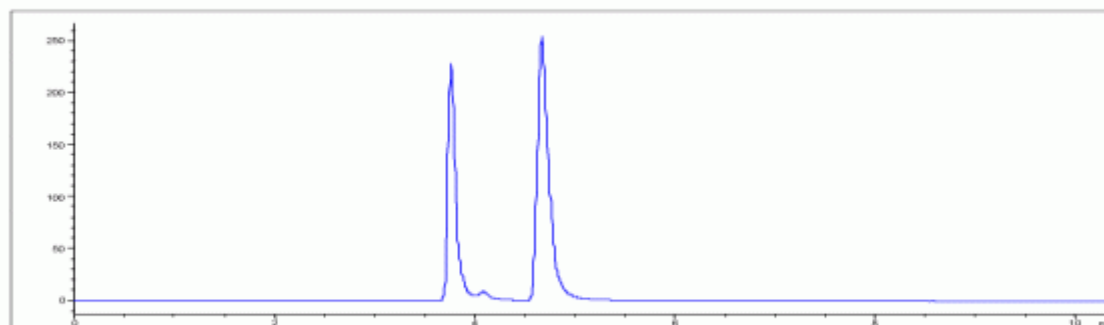

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	methanol / diethylamine = 100 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 260 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.7	0.24	1 - 2 2.26	5.52
2	4.6	0.54		